

USER MANUAL



MPW-260 MPW-260R

Read before use!

Serial number of centrifuges:

For centrifuges with serial no (SN):

MPW-260 10260045624 – 10260046524

MPW-260R 10260R171824 – ...

This manual was prepared with special care. MPW MED. INSTRUMENTS may change the manual at any time and without notice because of improvements, typographical errors or improvements to facilities. All rights reserved. No part of this User Manual may be modified, distributed, published or reproduced without the prior permission of MPW MED. INSTRUMENTS.

You can find the current version of the user manual on our website under:











mpw.pl/en – DOWNLOAD section (one can choose demanded language version of website).

Content


1. Symbols used in the manual and on the device	5
1.1 Markings on the device	5
2. Application	6
3. Technical specification	7
3.1 Environmental conditions	7
4. Installation	8
4.1 Content of the package	8
4.2 Location selection	8
4.3 Preparation for installation	9
4.4 Centrifuge installation	9
4.5 Starting the centrifuge	9
4.6 Opening and closing the cover	10
4.7 Current protection	10
5. Safety notes	10
5.1 General remarks	10
5.2 Placing the rotor and accessories in the centrifuge	11
5.3 Filling tubes	11
5.4 Filling the rotor	12
5.4.1 Angular rotors	12
5.4.2 Horizontal rotors	12
5.5 Safety hints	14
5.6 Operating conditions	14
5.7 Equipment life	15
5.8 Work safety	15
5.9 Unbalance	16
5.10 Emergency stop	16
5.11 Residual risk	16
5.12 Obligation to report a serious device incident	16
6. Product description	17
6.1 Product Design and Appearance	17
6.2 Name plate	18
6.3 Control device	18
6.4 Setting parameters	18
6.5 Safety features	18
6.6 Increase in temperature (MPW-260 only)	19
7. Centrifuging	20
7.1 Control panel	20
7.2 Display	20
7.3 Setting up RPM, RCF, time, temperature	22
7.4 User's programs	24
7.5 Creator of acceleration and deceleration curves	26
7.5.1 Acceleration characteristic, creation of episode 1	26
7.5.2 Adding and editing sections - acceleration	27
7.5.3 Acceleration graph	28
7.5.4 Deceleration characteristic, creation of episode 1	28
7.5.5 Adding and editing sections - deceleration	29
7.5.6 Deceleration graph	29
7.5.7 Deleting sections	30
7.6 Programs with user characteristics	30
7.7 Rotor and bucket choosing	30
7.8 SHORT mode	31
7.9 Finishing the centrifuging	31
7.10 Temporarily disabled functions	32
8. Temperature control	32



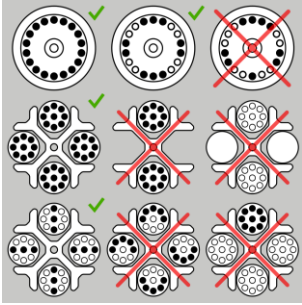


8.1	Initial cooling during centrifuging –FAST COOL	32
8.2	Initial cooling without centrifuging – THERMAL CHAMBER	33
8.3	Cooling in “START DELAY – OF TEMPERATURE” mode.....	33
8.4	Cooling in „SHORT” mode	33
8.5	Cooling notes.....	33
9.	Parameters of centrifugation	34
9.1	Acceleration/deceleration – changing characteristics	34
9.2	Radius	35
9.3	Sample density	35
9.4	Temperature offset	36
9.5	Thermal chamber	36
9.6	Automatic lid opening	37
9.7	Start delay - of time	37
9.8	Start delay – of temperature	38
10.	Menu	39
10.1	Screen saver	39
10.2	Visual alarm	40
10.3	Types of main screen.....	40
10.3.1	Switching the normal display to simplified display	40
10.3.2	Switching the simplified screen to normal display	41
10.4	Rotating runtime	42
10.5	Buzzer	42
10.6	Language	43
10.7	Other	43
10.8	Password	44
10.9	Last 10 cycles	45
10.10	Work time	46
10.11	Rotor runtime	46
10.12	Contact us	46
10.13	Diagnostics.....	47
10.14	Factory settings	47
11.	Maintenance	47
11.1	Cleaning of the centrifuge	47
11.2	Maintenance of centrifuge elements	48
11.3	Sterilization.....	48
11.3.1	Autoclaving	49
11.4	Chemical resistance	50
12.	Troubleshooting	51
12.1	Messages	52
12.2	Emergency cover release.....	53
13.	Guarantee	53
14.	Transport and storage	53
15.	Disposal	54
16.	List of changes in the manual	54
17.	Manufacturer’s info	55
	Distributor’s info	55
18.	ANNEXES	55
A.	Additional accessories	
B.	Declaration of conformity (CE, ROHS 2)	
C.	Declaration of decontamination (repair / return)	
D.	Nomogram RPM / RCF	

1. Symbols used in the manual and on the device

Symbol	Explanation
	WARNING! Warning of potential injury or health risk
	WARNING! Risk of crushing injury
	DANGER! Risk of electric shock with potential for severe injury or death as a consequence
	DANGER! Biohazard with potential for risk to health or death as a consequence
	DANGER! Risk of explosion with potential for severe injury or death as a consequence
	Symbol identifying a medical device for in vitro diagnostic use
	CE mark
	Symbol informing about the method of disposal
	Please read the instruction manual before you start working with the device
	Manufacturer's data

1.1 Markings on the device

Symbol	Explanation	Location
	Information about the direction of rotation of the rotor	Under the centrifuge lid

	<p>Information on where and how to use the emergency lid opening mechanism</p>	<p>On the side of the centrifuge next to the emergency opening of the lid</p>
	<p>Reminder for proper rotor maintenance</p>	<p>Under the centrifuge lid</p>
	<p>Information about correct and incorrect filling of rotors</p>	<p>Under the centrifuge lid</p>
 <p>Uwaga! Przed awaryjnym otwarciem pokrywy, wyłączyć urządzenie i odłączyć kabel zasilający. Odczekać 10 min i/lub zaglądnąć przez wziernik, upewnić się, że wirnik nie obraca się, a następnie otworzyć pokrywę.</p> <p>Attention! Before emergency opening the cover, switch off the mains power switch and disconnect the power cord. Wait 10 min and/or looking through the sight glass, make sure that the rotor is not rotating.</p>	<p>Information about the place of danger</p>	<p>On the side of the centrifuge next to the emergency opening of the lid</p>
 <p>CAUTION! UWAGA!</p> <p>Tighten the rotor fixing screw with the provided key. Dokrećić śrubę mocującą wirnik za pomocą dostarczonego klucza.</p>	<p>Information reminding about the proper tightening of the rotor screw</p>	<p>Under the centrifuge lid</p>

2. Application

- The **MPW-260, 260R** centrifuges are a family of table-top non-automatic laboratory centrifuges (MPW-260 – ventilated centrifuge, MPW-260R – refrigerated centrifuge).
- The devices are intended for In Vitro Diagnostics (IVD). This means that it is an in vitro diagnostic medical device - in accordance with the Regulation 2017/746 of the European Parliament and of the Council (EU) of 5 April 2017 on in vitro diagnostic medical devices and repealing Directive 98/79/EC and Commission Decision 2010 /227/EU.
- Centrifuges are used to separate aqueous solutions and suspensions of samples with a density not higher than **1.2g/cm³** taken from human, animal and plant organisms into components of different densities under the influence of centrifugal force, in order to provide information about their biological state and to other analytical work.
- The design of the centrifuge ensures ease of use, safe operation and a wide range of applications in medical, biochemical, and other analysis laboratories.
- The centrifuge is not biotight, therefore, when centrifuging preparations that require biotightness, containers and rotors with a biotightness certificate should be used.

3. Technical specification

manufacturer	"MPW MED. INSTRUMENTS" SPÓŁDZIELNIA PRACY, Boremłowska 46 Street, 04-347 Warsaw										
type	MPW - 260					MPW - 260R					
cat. number (REF)	10260/2-56	10260/ 1-56/100	10260/ 1-56/110	10260/1-56	10260/ 1-56/127	10260R/2-5	10260R/2-6	10260R/ 1-6/100	10260R/ 1-6/110	10260R/1-6	10260R/ 1-6/127
mains voltage (L1+N+PE)	230V	100V	110V	120V	127V	230V	230V	100V	110V	120V	127V
	±10%	±5%				±10%		±5%			
mains frequency,	50/60Hz	50Hz/60Hz				50Hz	60Hz	60Hz			
current protection [A]	T 4A	T 6,3A				T 10A					
cooling medium						R452A					
power consumption	230W					530W					
	MPW - 260					MPW - 260R					
capacity (max.)	500 ml										
speed – RPM	90 ÷ 18000 rpm (step 1 rpm)										
force – RCF	24270 x g (step 1 x g)										
kinetic energy (max.)	11000 J										
running time	00:00:01 ÷ 99:59:59 – [hours, min., sec] (step 1s)										
time counting	since start button is pressed / since preselected speed is reached										
short-time operation mode – SHORT	yes										
continuous operation mode – HOLD	yes										
menu languages	Polish, English, German, Spanish, Italian, Portuguese, Russian, Swedish, French, Czech										
user programs	100										
adjustable temperature						-20 ÷ 40°C* (step 1°C)					
guaranteed temperature with max. rotor speed						≤4°C					
cooling/heating without centrifuging	no					yes/ no					
cooling/heating with centrifuging	no					yes /no					
acceleration (ACEL)	10 linear curves										
deceleration (DECEL)	10 linear curves										
programmable non-linear curves:											
acceleration	10										
deceleration	10										
USB communication	no										
Electromagnetic compatibility	according to EN 61326-2-6:2006										
degree of protection (according to PN-EN 60034-5:2021-01)	IP20										
height (H)	315 mm					315 mm					
width (W)	365 mm					365 mm					
depth (D)	495 mm					660 mm					
height with open cover (H _{oc})	620 mm					620 mm					
noise level	<60 dB										
weight 230V	approx. 23,1 kg					approx. 42,2 kg					
weight 120V	approx. 25 kg					approx. 44,3 kg					

*time and possibility of obtaining a set temperature is dependent on multiple factors, including rotor type, established RPM, ambient temperature, accuracy ±1°C, appropriate for place of temperature sensor

3.1 Environmental conditions

- The device may only be used indoors.
- The permissible ambient temperature is 2°C to 40°C.
- Maximum allowed relative humidity 80% at temperature up to 31°C decreasing linearly to 50% relative humidity at 40°C.
- The mains voltage fluctuations must not exceed ± 10% of the nominal voltage.
- Maximum altitude 2,000 m above sea level.


- Overvoltage category II.
- Pollution degree 2.

4. Installation


4.1 Content of the package



name	pcs.	cat no.
centrifuge MPW-260/R (type and supply version dependent)	1	10260/2-56; 10260/1-56; 10260/1-56/100; 10260/1-56/110; 10260/1-56/127; 10260R/2-5; 10260R/2-6; 10260R/1-6; 10260R/1-6/100; 10260R/1-6/110; 10260R/1-6/127
complete clamp	1	17142
spanner for the rotor	1	17099T
spanner for emergency opening of the cover	1	18640
power cord – 230V / 120V	1	17866/17867
fuse WTA T 4A – 230V / WTA T 6,3A – 120V (for MPW-260)	2	17861 17862
fuse WTA T 10A – 230V / 120V (for MPW-260R)	2	17863
vaseline 20ml	1	17201
user manual	1	See page 1.

4.2 Location selection

	<p>WARNING! Risk of damage to the device.</p> <ul style="list-style-type: none"> ▪ The table intended for the centrifuge operation should be adapted to the weight of the device, clean, stable and free from vibrations, and have a flat, levelled top. ▪ In accordance with the EN 61010-2-020 standard, leave a safety distance of 30 cm from the operating device. Keep a distance from walls and other devices. Do not place any objects in this area. ▪ The centrifuge should be positioned so that access to the mains switch is not difficult. ▪ Do not use the device near strong, unshielded, high-frequency electromagnetic sources as they may interfere with its proper operation. ▪ Do not install the centrifuge near heat sources (e.g. radiators). ▪ Avoid direct sunlight. ▪ Ensure adequate ventilation of the room.
---	---

4.3 Preparation for installation

	<p>WARNING! Risk of injury or damage to the device.</p> <ul style="list-style-type: none">▪ After changing the storage location of the device (from cold to warm), wait until the device warms up to ambient temperature to avoid damage to electronic components due to condensation.▪ It is important to allow enough time for the device to dry before restarting it (min. 4 hours).▪ Lifting and carrying the device may result in injuries due to its heavy weight.▪ The centrifuge should be lifted and transported with a sufficient number of people (min. 2). Use a transport aid to move the centrifuge.▪ Lift the device from below, near its feet.
---	---

 	<p>WARNING! Risk of electric shock or fire.</p> <ul style="list-style-type: none">▪ The centrifuge may only be operated in a building that complies with applicable national regulations and standards. In particular, it must be ensured that power supply circuits located upstream of the device's internal protection are not loaded in an unauthorized manner. This can be ensured by using additional interrupters or other suitable fuse elements in the building installation.▪ The voltage and frequency of the power source must comply with the requirements specified on the device nameplate.▪ The power socket must be earthed with a protective conductor (PE).▪ During operation, there must be easy access to the power switch and the device that cuts off the electrical network (e.g. residual current device).▪ Only the power cord included with the centrifuge can be used.▪ Before turning on the device, make sure it is properly connected to the power source.
--	--

- 1) Open the package.
- 2) Remove the box containing the accessories.
- 3) Remove the centrifuge from the box and remove the foil (keep the packaging and packing material for service shipment).
- 4) Place the device on a suitable laboratory table.

4.4 Centrifuge installation


- 1) Check whether the mains voltage and frequency meet the requirements given on the nameplate of the device.
- 2) Connect the power cord to the centrifuge power socket (on the rear wall of the centrifuge) and to the power source.


4.5 Starting the centrifuge


- 1) Wait at least 4 hours for the unit to reach ambient temperature to avoid compressor failure or damage to electronic components due to condensation.
- 2) Turn on the centrifuge power using the mains switch located on the side wall of the device.
- 3) Open the cover according to the section **Opening and closing the cover**.
- 4) Install the rotor according to the section **Placing the rotor and accessories in the centrifuge**.

- 5) Set centrifugation parameters according to the sections **Centrifuging** and **Parameters of centrifugation**.


4.6 Opening and closing the cover

	<p>ATTENTION !</p> <ul style="list-style-type: none"> ▪ The cover can only be opened when the centrifuge is at rest (the rotor is not rotating). ▪ Centrifugation can only be started with the lid closed.
---	---

	<p>WARNING! Risk of injury.</p> <ul style="list-style-type: none"> ▪ Do not put your hands between the cover and the housing when closing the centrifuge cover.
---	---


- 1) Press the **COVER** button  to open the cover.
- 2) To close the lid, press it down with both hands until the lock engages.

4.7 Current protection


	The centrifuge is equipped with thermal current protection. Fuse is situated in the plug-in socket unit at back wall of the centrifuge.
---	---

5. Safety notes


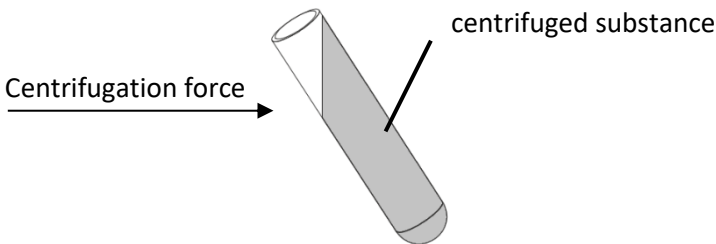
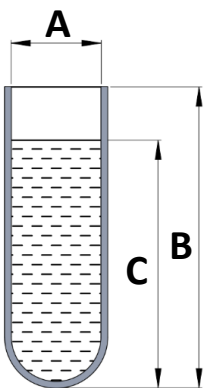
5.1 General remarks

	<ul style="list-style-type: none"> ▪ The laboratory centrifuge may be operated only by qualified laboratory personnel, after reading the operating manual. ▪ The operating instructions are part of the product. ▪ The operating manual should always be kept in the vicinity of the centrifuge. ▪ The centrifuge cannot be operated contrary to its purpose. ▪ If the centrifuge is used in a manner inconsistent with the manufacturer's guidelines, the safety of its use may be impaired. ▪ For centrifugation in the centrifuge, only containers and inserts provided in the list of equipment and centrifuge tubes, the diameter, length and strength of which are appropriate, should be used. The use of test tubes not included in the list should be agreed with MPW MED. INSTRUMENTS or its authorized representatives. ▪ Pay attention to the quality and appropriate thickness of the glass test tubes walls. Glass tubes should be centrifuge tubes, and their use in the centrifuge should be made dependent on the following guidelines: <table border="1" style="width: 100%; border-collapse: collapse;"> <thead> <tr> <th style="text-align: left;">glass tubes</th> <th style="text-align: center;">max RCF in angular rotors</th> <th style="text-align: center;">max RCF in horizontal rotors</th> </tr> </thead> <tbody> <tr> <td>5-10 ml</td> <td style="text-align: center;">3000 x g</td> <td style="text-align: center;">4000 x g</td> </tr> <tr> <td>30-100 ml</td> <td style="text-align: center;">spinning not allowed</td> <td style="text-align: center;">4000 x g</td> </tr> </tbody> </table> <ul style="list-style-type: none"> ▪ Weighing the filled test tubes into the rotor is recommended. When centrifuging in horizontal rotors, it is recommended to weigh the filled containers / hangers. This will allow to minimize the differences in mass between them, and as a result to avoid the negative impact of vibrations on the engine suspension and to reduce noise levels during the operation of the centrifuge. 	glass tubes	max RCF in angular rotors	max RCF in horizontal rotors	5-10 ml	3000 x g	4000 x g	30-100 ml	spinning not allowed	4000 x g
glass tubes	max RCF in angular rotors	max RCF in horizontal rotors								
5-10 ml	3000 x g	4000 x g								
30-100 ml	spinning not allowed	4000 x g								

5.2 Placing the rotor and accessories in the centrifuge



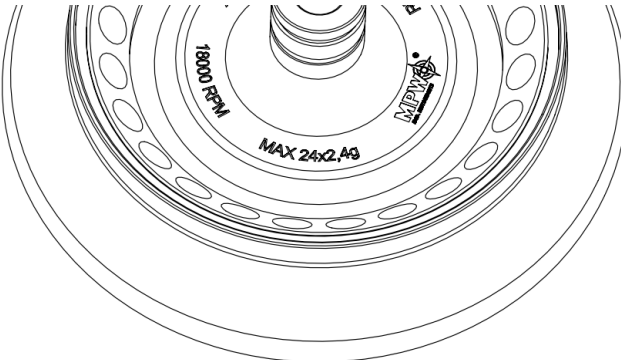
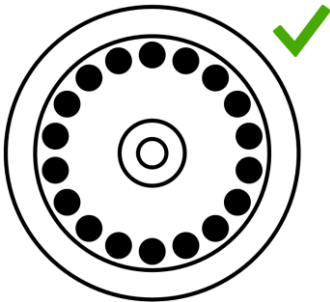
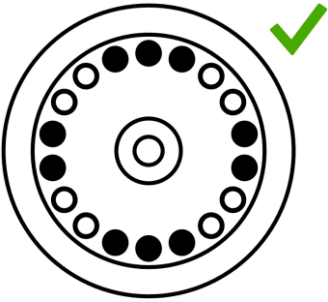
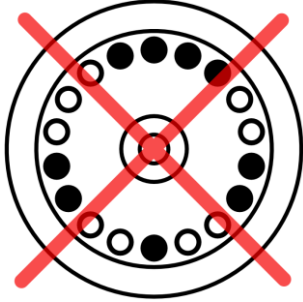
	<ul style="list-style-type: none"> ▪ Connect the centrifuge to the power supply (mains socket at the back of the centrifuge). ▪ Turn on the centrifuge (switch on the side of the centrifuge). ▪ Open the cover of the centrifuge by pressing the COVER key. Before installing the rotor, check that the rotating chamber is free from any contamination. If there is dust, glass splinters, liquid residues, etc., remove them. ▪ The rotor can fall if not handled properly, therefore it should always be handled and placed in the centrifuge using both hands. ▪ Place the rotor on the motor axis by sliding it onto the cone as far as it will go (keeping the coaxiality between the rotor and the motor axis). ▪ Screw the screw fixing the rotor into the motor axis (clockwise), and then tighten it firmly with the rotor key. ▪ Fill the rotor with containers / hangers / test tubes according to recommendations in section Filling the rotor. ▪ In order to replace the rotor, first remove the tubes and containers from it, unscrew the screw fixing the rotor with the enclosed rotor key, counterclockwise, then using both hands, grab the rotor on opposite sides and remove it from the motor axis. ▪ Install another rotor as described above instructions.
---	---

5.3 Filling tubes

	<ul style="list-style-type: none"> ▪ Fill test tubes outside the centrifuge. <div style="text-align: center; margin: 10px 0;">  </div> <ul style="list-style-type: none"> ▪ If the manufacturer of the test tube has not specified the maximum level, fill the test tubes so that the centrifuged substance does not run out of the vessel during centrifugation. To do this, use the formula below: <div style="display: flex; align-items: center; justify-content: center; margin: 10px 0;"> <div style="text-align: center;">  </div> <div style="margin-left: 20px;"> $C < B - \frac{A}{2}$ </div> </div> <div style="margin-top: 10px;"> <p>A – internal tube diameter B – tube height C – max liquid level</p> </div>
---	--

5.4 Filling the rotor

5.4.1 Angular rotors

	<p>CAUTION!</p> <ul style="list-style-type: none"> Angle rotors must be used with a suitable cover which must be screwed securely onto the rotor. The rotor and the cover are marked with the same catalog number (REF) to eliminate the risk of incorrect selection when you have several types of rotors. 	
	<ul style="list-style-type: none"> Check that the impeller is seated correctly and firmly bolted to the motor shaft. Do not exceed the maximum rotor load (information is provided on the rotor). <p style="text-align: center;">An example of the marking on the angular rotor:</p>  <p style="text-align: center;">MAX. 24x2,4g - means the possibility of placing 24 test tubes in the rotor, each weighing 2.4 g.</p> <ul style="list-style-type: none"> To ensure symmetrical loading, insert test tubes of the same type and weight in pairs into opposite openings of the rotor. If reduction inserts are used, they should also be placed in the holes opposite to each other in pairs of the same type. 	
<p>Examples of correct and incorrect arrangement of test tubes in the rotor:</p>		
		

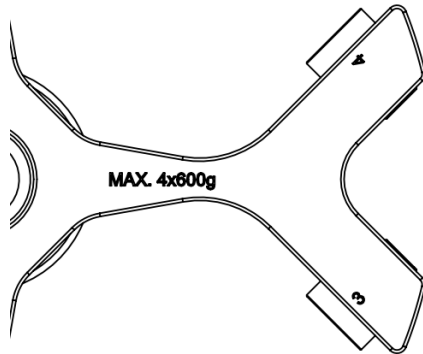
5.4.2 Horizontal rotors

- Check that the impeller is seated correctly and firmly bolted to the motor shaft.
- Make sure that the rotor pins and grooves of the containers / hangers are clean, and then it is necessary to lubricate them with the technical petroleum jelly supplied with the device (catalog number 17201).
- Place the containers / hangers in the rotor.
- Horizontal rotors must be filled with a set of containers / hangers.
- Observe the limitations for the permissible centrifugal mass stated on the rotor and container. If the marking appears on the rotor, it refers to the mass of the substance to be centrifuged,

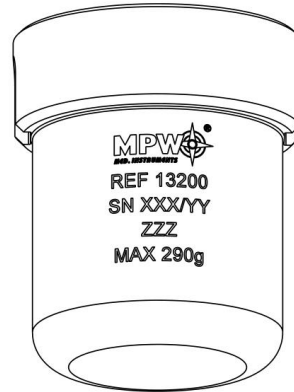
and if on the container it refers to the mass of the contents of the container, i.e. insert, test tube and the substance contained in it.

Examples of markings on horizontal rotors and containers:

Marking on the rotor



Marking on the container

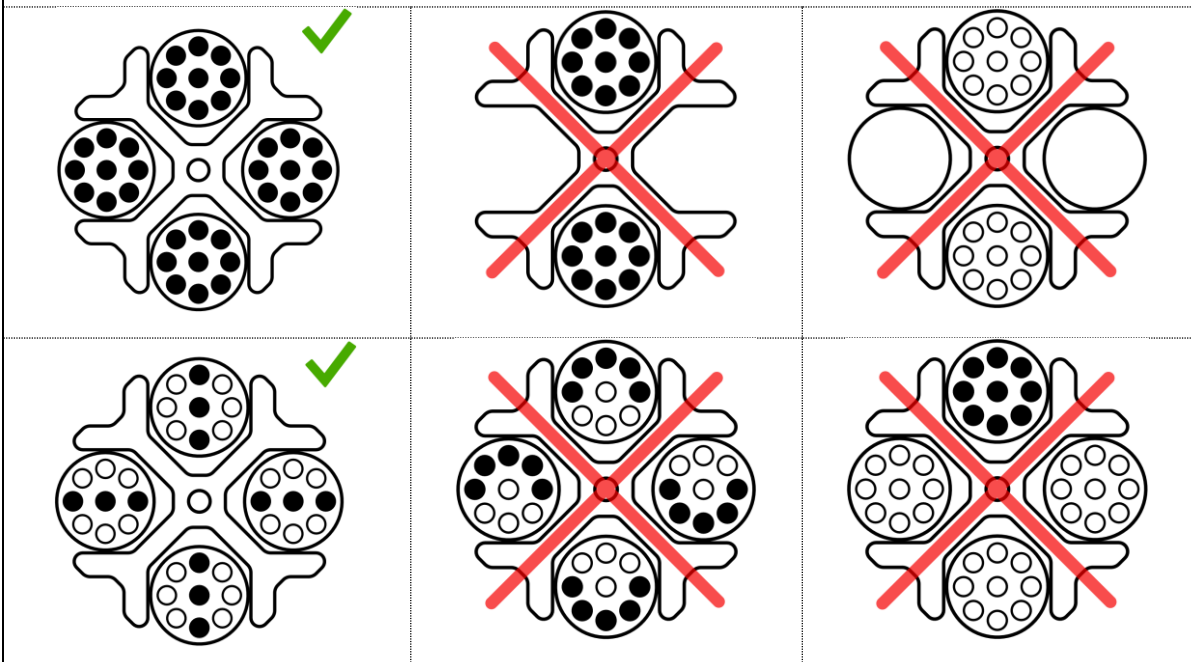


MAX. 4x600g – permissible weight of the contents of the test tubes placed in each of the 4 containers





MAX. 290g – maximum weight of the contents of the container

- In order to ensure symmetrical and even rotor load, try to fill opposite seats with containers / hangers of the same type and weight.
- Tubes should be placed symmetrically facing each other.
- Before starting the centrifugation, check that all containers / hangers are properly placed in the rotor and can swing freely.
- For this purpose, empty test tubes should be placed in containers. Manually tilt the containers to the horizontal position and check that there are no collisions between the tubes, containers / hangers and the rotor.



Examples of correct and incorrect arrangement of test tubes in the rotor:




5.5 Safety hints

	<p>ROTOR MAINTENANCE</p> <ul style="list-style-type: none"> ▪ In order to increase the durability of gaskets, threaded places, rotor pins, undercuts for pins in containers, they must be cleaned, and then it is necessary to lubricate them with the technical petroleum jelly supplied with the device (catalog number 17201). ▪ Use only accessories that are in good technical condition.
	<p>HU EQUIPMENT MAINTENANCE</p> <ul style="list-style-type: none"> ▪ Make sure the sealing rings (rubber) are lightly coated with grease to maintain tightness. Use high vacuum silicone grease, eg type "C" by LUBRINA.
	<p>HAZARDOUS MATERIALS</p> <ul style="list-style-type: none"> ▪ Infectious materials should be centrifuged only in containers / rotors with covers. ▪ It is not allowed to centrifuge toxic or infectious materials if the rotor or test tube seal is damaged. ▪ Appropriate disinfection procedures should always be carried out, if hazardous substances have contaminated the centrifuge or its accessories.
	<p>EXPLOSIVE, FLAMMABLE MATERIALS</p> <ul style="list-style-type: none"> ▪ It is not allowed to centrifuge explosive and inflammable materials. ▪ Do not centrifuge substances that could create a potentially explosive atmosphere as a result of the high energy supply during centrifugation. ▪ The centrifuge must not be used in an explosive atmosphere. ▪ It is not allowed to centrifuge materials that may generate flammable or explosive mixtures when exposed to air.

5.6 Operating conditions

	<p>GENERAL REMARKS</p> <ul style="list-style-type: none"> ▪ Only original equipment of centrifuges and spare parts should be used. ▪ In case of a malfunction of the centrifuge, the MPW MED factory service should be used. INSTRUMENTS or its authorized representatives. ▪ It is not allowed to start the centrifuge if it is not installed correctly or the rotor and accessories are not properly mounted. ▪ The centrifuge must not be transported with the rotor installed on the motor shaft. ▪ Fill the rotor equipment to the same weight in order to prevent unbalance of the centrifuge (point Filling the rotor).
	<p>START-UP</p> <ul style="list-style-type: none"> ▪ Before switching on the device, carefully read all sections of this manual in order to ensure the correct operation of the device and to avoid damage to the device or its accessories.

	<p>CENTRIFUGAL SUBSTANCES</p> <ul style="list-style-type: none"> ▪ Rotors are designed for centrifuging liquids with an average density of 1.2 g / cm³ or less. This applies to centrifugation at maximum speed. If liquids with a higher density are to be used, be sure to enter the density value in the PARAM / DENSITY tab in order to reduce the available spin speed.
---	--




5.7 Equipment life

	<ul style="list-style-type: none"> ▪ Each spin cycle in which the rotor has accelerated and decelerated is considered a duty cycle, independent of speed and duration. ▪ Do not use the equipment after the allowable number of cycles or after the maximum service life has passed, whichever comes first (the service life is 15000 cycles or 5 years).
--	--

5.8 Work safety


The centrifuge should be inspected by an authorized service at least once a year (after the warranty period). Special circumstances, e.g., corrosive environment, may be the reason for more frequent checks. Tests should end with issuing a validation protocol, which specifies checking the technical condition of a laboratory centrifuge.

It is recommended to create a document that records all repairs and inspections. This document should be kept in the place where the centrifuge is used.

	<p>CONTROLS CONDUCTED BY THE OPERATOR</p> <ul style="list-style-type: none"> ▪ The operator must pay attention to the fact that the parts of the centrifuge, important from the safety point of view, are not damaged. This remark applies to: <ul style="list-style-type: none"> ▪ Centrifuge accessories, especially structural changes, corrosion, initial cracks, abrasion of metal parts. ▪ Bolted connections. ▪ Inspection of rotor and container seals, if any. Particular attention should be paid to rubber elements (seals). In the event of any damage or visible structural changes, they should be immediately replaced with new ones. ▪ Control of the performance of annual post-warranty inspections of the technical condition of the centrifuge. ▪ During centrifugation, it is not allowed to lift, shift the centrifuge or rest on it. ▪ During centrifugation one must not stay in the safety zone, i.e., 30 cm distance around the centrifuge, nor leave any objects, e.g., glass vessels, inside this zone. ▪ It is not allowed to put any objects on the centrifuge.
	<p>OPENING THE COVER DURING SPINNING</p> <ul style="list-style-type: none"> ▪ It is not allowed to use the emergency cover opening during centrifuging, because it may result in loss of health or life.
	<p>HANDLING OF ROTORS</p> <ul style="list-style-type: none"> ▪ It is not allowed to use accessories (rotors, lids, containers, hangers and round carriers) with signs of corrosion or other mechanical damage.

	<ul style="list-style-type: none"> ▪ It is not allowed to centrifuge substances of high corrosive aggressiveness, which may damage the materials and reduce the mechanical properties of rotors, buckets and round carriers. ▪ It is not allowed to centrifuge rotors with removed or loose covers.
--	---

5.9 *Unbalance*

	<p>Unbalance causes noise, vibration during operation and has a negative effect on the driveline (engine and suspension). The more precisely the process of balancing the feed to the rotor is carried out, the smoother the centrifuge will run and the longer the useful life of the drive system will be. Moreover, thanks to the correct balancing, an excellent level of separation of the centrifuged substance is achieved since the separated components will not be picked up again by vibrations.</p>
---	---

The centrifuge is equipped with a rotor imbalance sensor. In the event of its activation, the centrifugation process is stopped by quick braking and an error message is displayed. Erasing the error message is possible by pressing one of the following buttons: **BACK, STOP, COVER, SET** and **▲ ▼ ◀▶**.

Make sure that the rotor has been properly loaded - places in the rotor must be equipped with identically filled containers, inserts and test tubes so as to obtain the best possible weight balance (see chapter Filling the rotor). If necessary, correct the load distribution and / or, in the case of horizontal rotors, clean and lubricate the rotor pins, then restart the spin.

5.10 *Emergency stop*

At any time during centrifugation, it is possible to interrupt the process and stop the centrifugation with the fastest rotor characteristics. This is done by pressing the stop button twice (**2x STOP**).

Pressing the **STOP** key once will stop the spinning with the braking characteristics set in the program. The message about interrupted centrifuging can be canceled with the following buttons: **BACK, STOP, COVER, SET** and **▲ ▼ ◀▶**.

5.11 *Residual risk*

The centrifuge is built according to the state-of-the-art and the recognized safety regulations. Nevertheless, still remain some level of residual risk due to improper operation and malfunctions. It is possible to decrease residual risk by strictly applying user manual conditions and correcting malfunction which could threaten safety, immediately.

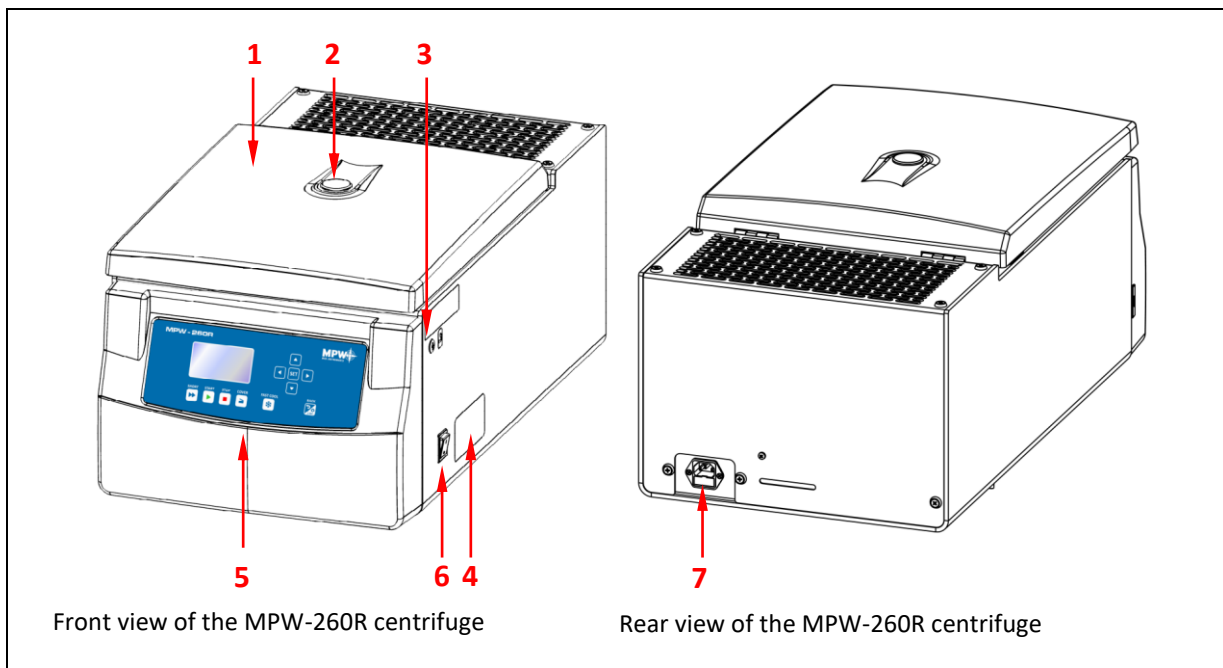
5.12 *Obligation to report a serious device incident*

Any serious incident related to the device should be reported to the manufacturer and the competent authority of the Member State where the user or patient resides.

6. Product description

6.1 Product Design and Appearance

A new generation of MPW MED laboratory centrifuges. INSTRUMENTS is equipped with modern microprocessor controllers, very durable and quiet brushless asynchronous motors and equipment that meets modern user requirements. The centrifuge has a rigid self-supporting structure. The housing is made of lacquered aluminum sheet, the back is made of steel sheet. The front part and the cover are made of ABS plastic. The cover is mounted on steel hinges, and from the front it is secured against opening during spinning with an electromagnetic lock. The spinning chamber is made of stainless steel.



- 1 Centrifuge cover
- 2 Inspection glass (to control the rest condition of the rotor)
- 3 Place for emergency opening of the cover
- 4 Name plate
- 5 Control panel (display and control of the centrifuge)
- 6 Main's switch
- 7 Centrifuge power socket (with fuse's socket)

6.2 Name plate

Data regarding the device should be read from the nameplate located on the side wall of the centrifuge next to the power switch (the image below is an example).

1	Model	MPW-260R
2	REF	10260R/2-5
3	n	18000 rpm
4	AC	230 V
5	P	530 W
6	Kinetic energy	11000 J
7	SN	10260R000124
8	Approval marks and symbols	CE, IVD, i, warning, crossed-out trash
9	Rated frequency	50 Hz
10	Current protection	T 10 A
11	Manufacturer's logo	MPW
12	Manufacturer's information	MPW MED. INSTRUMENTS SPÓŁDZIELNIA PRACY ul. Boremłowska 46, 04-347 Warszawa, POLAND
13	Information about the refrigerant	R452A / 0,175 kg / 0,38 t eq CO ₂ / GWP 2141 / 27 bar fluorinated greenhouse gas in hermetically sealed refrigeration system
14	Density of centrifuged substance	1,2 g/ml
15	UDI code carrier – data matrix (AIDC) with human-readable presentation (HRI)	UDI (01)5905386360272 (11)20240108 (21)10260R000124
16	Date of production	2024-01-08

1	Centrifuge model
2	Catalog number
3	Maximum speed
4	Rated voltage
5	Maximum rated power
6	Kinetic energy
7	Serial number
8	Approval marks and symbols (explained in chapter 1)
9	Rated frequency
10	Current protection
11	Manufacturer's logo
12	Manufacturer's information
13	Information about the refrigerant (refrigerated centrifuges only)
14	Density of centrifuged substance
15	UDI code carrier – data matrix (AIDC) with human-readable presentation (HRI)
16	Date of production

6.3 Control device

The microprocessor control unit of the centrifuge ensures broad possibilities of providing, realization and reading of work parameters.

6.4 Setting parameters

Data setting and read-out system forms hermetically closed keyboard with distinctly accessible operation points. Easily readable displays signaling individual performed operations facilitate operator's programming and recording of parameters and condition of the centrifuge.

6.5 Safety features

Cover lock

The centrifuge can be started only with properly closed cover. While the cover can be opened only after stopping the rotor. In case of emergency opening of the cover during operation, the centrifuge drive will be immediately switched-off and the rotor will brake till complete stopping.

Unbalance detecting




When loads of opposite buckets or carriers in rotors are unbalanced, the drive will be switched-off during acceleration or operation of the centrifuge – and the error message will be displayed.

Rotor verification and checking compatibility with loaded program

Directly after starting centrifuging, a unit verifies the type of the rotor applied and in the case of its incompatibility with the type indicated in the application or absence of the rotor, the spinning process shall be stopped with simultaneous displaying the error message. The conformity of the type of the

rotor is signaled with a single audible signal. In case auto identification (see “9.8 Other”) option is checked, proper rotor will be automatically chosen, without user engagement.

Rest state inspection

Opening of the centrifuge’s cover by **COVER** button is possible only when the rotor is in the state of rest. Check if the symbol  is visible on the screen. Use inspection glass in cover for be sure if rotor is in the rest state. When the rotor is being stopped, braking symbol  or  (see chapter “Display”) is visible and goes off when it is stopped. Emergency cover opening during rotor running is prohibited.

Checking of excessive temperature (MPW-260R)

If temperature in rotation chamber exceeds 50°C caused by, for example, malfunction of cooling system, drive will be switched off and error message will be displayed. The reboot is only possible after chilling device.

6.6 Increase in temperature (MPW-260 only)

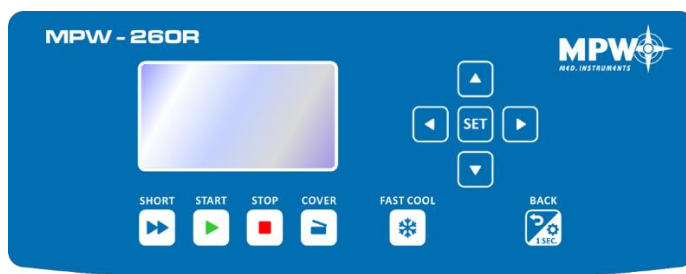
In uncooled centrifuges, the temperature in the rotor chamber, rotor and sample can increase to above 40°C, based on the run time, g-force (RCF)/speed and ambient temperature.

7. Centrifuging

Power switching ON/OFF is carried out with master switch situated on the right-side wall of the centrifuge. All settings on the centrifuge are done by means of the control panel.

7.1 Control panel

The control panel placed on the front casing serves the purpose of controlling centrifuge operation.



Control panel MPW-260R

	SHORT ¹	short-time centrifuging
	START	start centrifugation run
	STOP ²	end centrifugation run
	COVER	cover opening
	FAST COOL	start fast cooling mode (MPW-260R only)
	BACK/ OPTIONS	exit the current menu / enter to submenu of options (keep held down within 1 s.)
	UP	navigation in menu / increasing values
	DOWN	navigation in menu / decreasing values
	LEFT	navigation in menu
	RIGHT	navigation in menu
SET	SET	changing parameters / confirming changes

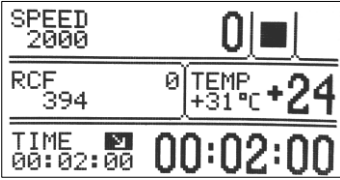
¹ the centrifuge is working as long as the key is pressed

² first-time pressing press – will make stopping centrifuging with acceleration characteristics set in the current program,
second-time pressing – will make stopping the centrifuging as fast as possible (quickest characteristic) (after stopping the rotor, the message can be cancelled by pressing any key except **SHORT**, **START** and **COVER** – if cover is open)


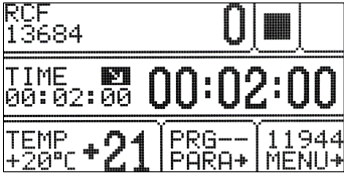
7.2 Display

The display is located in the center of the control panel. The main screen variants are presented below. In this chapter are shown exemplary screens for MPW-260R centrifuge (screens for MPW-260 – without cooling – do not include temperature field). Highlighting of field on display means it is selected and ready to set.

	<p>After switching on centrifuge, welcome screen appears. After disappearing the welcome screen, it is possible to setting up parameters.</p>
--	---

	<p>Simplified display mode is set as default, there is possible to switch to normal (see chapter <i>Types of main screen</i>).</p>
---	--








Normal display contains an expanded number of settings visible during operation.

RPM display mode	RCF display mode
	

- Detailed information on display modes is provided in chapter "Types of main screen".

SPEED	rotor speed	assigned/measured
RCF	relative centrifugal force	assigned/measured
TIME	centrifuging time	assigned/measured
TEMP	temperature	assigned/measured
PRG	program no.	
11199	rotor no.	
PARAM	parameters of the centrifuge	
MENU	configuration menu	

	changing values		
	user multi sections curve		
	density > 1,2 g/cm ³		
	centrifuging radius changed		
	counting time down (decreasing)		counting time up (increasing)
	cooling to assigned temperature		
	FAST COOL mode cooling		
	centrifuging		centrifuging (with automatic cover opening)
	rotor stopped / closed cover		rotor stopped / opened lid
	braking		fastest decelerating
	rotor identification		
	thermal chamber		
	temperature delay		

	time delay		
	drop-down list		
	temporarily disabled		
	locked		
	time counting (blinking)		
	disabled option		active option



7.3 Setting up RPM, RCF, time, temperature

On the main screen, it is possible to set:

rotating speed - RPM	SPEED
relative centrifugal force (multiple of g-force)	RCF
centrifuging time	TIME
centrifuging temperature	TEMP (only MPW-260R)

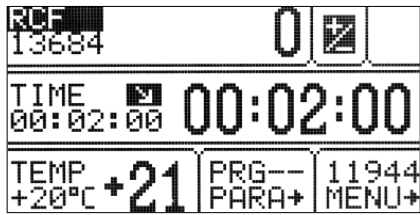
Exemplary change of **SPEED** setting:



- Press **SET** (to enter edit mode) –  appears.
- Via **▲▼◀▶** keys mark **SPEED** field (blinking).
- Press **SET** –  blinking.
- Via **◀▶** choose order of magnitude of changing value (blinking).
- With **▲▼** choose demanded value.
- Repeat above two steps for other orders of magnitude.
- Confirm settings by pressing **SET**.
- Press **BACK**.

- When RPM is changed, RCF is automatically corrected.

Exemplary change of **RCF** setting:

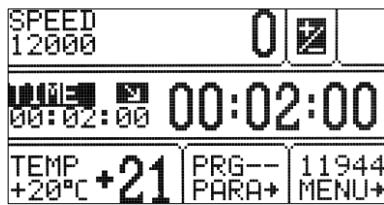


- Press **SET** (to enter edit mode) – appears.
- Via **▲▼◀▶** keys mark **RCF** field (blinking).
- Press **SET-** blinking.
- Via **◀▶** choose order of magnitude of changing value (blinking).
- With **▲▼** choose demanded value.
- Repeat above two steps for other orders of magnitude.
- Confirm settings by pressing **SET**.
- Press **BACK**.

- When RCF is changed, RPM is automatically corrected.

- When setting the speed value, setting "hundreds" or "thousands" resets the "units" and "tens".

Exemplary change of **TIME** setting:



- Press **SET** (to enter edit mode) - appears.
- Via **▲▼◀▶** keys mark **TIME** field (blinking).

00:02:00
[hh : mm : ss]

e.g.:
centrifuging time – 2 minutes 00 seconds

- Press **SET-** blinking.
- Via **◀▶** choose order of magnitude of changing value (blinking).
- With **▲▼** choose demanded value.
- Repeat above two steps for other orders of magnitude.
- Confirm settings by pressing **SET**.
- Exit edit mode by pressing **BACK**.

00:02:00

set value

02:00

current value (most significant digits)

HOLD mode – continuous run mode



- To run centrifuging in **HOLD** mode set **00:00:00** time.
- To end centrifuging in **HOLD** mode press **STOP**.

Exemplary change of **TEMP** setting (only MPW-260R):

SPEED 12000 0

TIME 00:02:00 00:02:00

TEMP +21 PRG-- 11944
+20°C PARA+ MENU+

- Press **SET** (to enter edit mode) – appears.
- Via **▲▼◀▶** keys mark **TEMP** field (blinking).
- Press **SET** key.
- With **▲▼** choose demanded value [-20°C÷40°C].
- Confirm settings by pressing **SET**.
- Press **BACK**.

7.4 User's programs

SPEED 12000 0

TIME 02:10:00 02:10:00

TEMP +21 PRG 2 11716
+15°C PARA+ MENU+

After switching centrifuge on, program that was used in previous session is being loaded. If any program was not used in previous session, centrifuge resume the last chosen parameters.

SPEED 12000 0

TIME 02:10:00 02:10:00

TEMP +21 PRG-- 11716
+15°C PARA+ MENU+

Program choosing:

Entering the program selection mode for the simplified display:

SP 2

RC .9

TI 00

PROG.-- 11199
PARAM.+
MENU+

DISPLAY MODE

- Press and hold by 1 second.
- Choose **PROG** with **▲▼**
- Press **SET**.

Entering the program selection mode for the normal display:

SPEED 12000 0

TIME 02:10:00 02:10:00


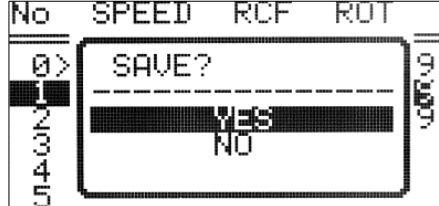
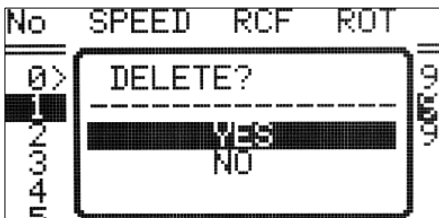
TEMP +21 PRG-- 11716
+15°C PARA+ MENU+

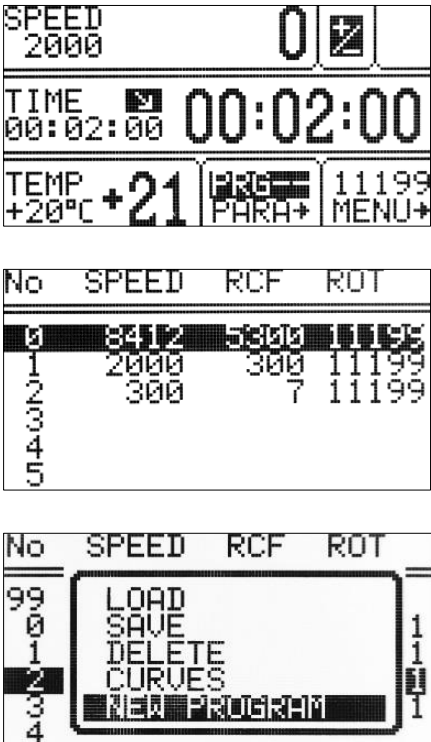
- Press **SET** key – appears.
- Via **▲▼◀▶** keys mark **PRG --** field (highlighted)
- Press **SET** key – list of programs is visible.

Program selection mode tab:

	SPEED	RCF	ROT
0	8412	5300	11199
1	2000	300	11199
2	300	7	11199
3			
4			
5			

- Via **▲▼** choose demanded program.
- Confirm with **SET** key.

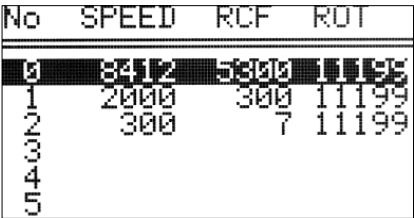


	<p>LOAD, SAVE, DELETE, CURVES, NEW PROGRAM refer chosen program which is marked by ▶.</p> <p>▶ – currently chosen program.</p>
	<ul style="list-style-type: none"> ▪ LOAD – load selected program
	<ul style="list-style-type: none"> ▪ SAVE – save settings as a program (confirm by selecting YES and pressing SET)
	<ul style="list-style-type: none"> ▪ DELETE – delete program (confirm by selecting YES and pressing SET)
	<ul style="list-style-type: none"> ▪ CURVES – creating characteristics
	<ul style="list-style-type: none"> ▪ NEW PROGRAM – creating new program

<p>NEW PROGRAM – enter to create new program mode (as below)</p>	
<p>Creating a new program:</p>	
	<ul style="list-style-type: none"> ▪ Press SET key. ▪ Via ▲▼◀▶ keys mark PRG -- field (blinking). ▪ Press SET key. List of programs is visible, choose demanded position (number of program). ▪ Press SET key- menu of program settings will appear. ▪ Choose NEW PROGRAM press SET and BACK, and then set demanded parameters of centrifuging (look chapter “6. Centrifuging”). ▪ In case you want to register new program, back to the PRG -- menu and save it as described before.



Changing parameters during centrifuging:

- There is a possibility to change parameters: **SPEED, RCF, TIME, TEMP** during centrifuging. Such modifications inactivate currently running program. When program was set, modification during run is represented by **PRG --** symbol (instead of the program number).

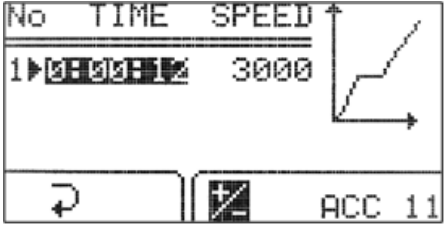

7.5 Creator of acceleration and deceleration curves

		PROG/CURVES
		<ul style="list-style-type: none"> With ▲▼ keys choose saved program for which you intend to create the acceleration or deceleration characteristics (marked with symbol ▶). Press SET. With ▲▼ keys choose CURVES. Press SET - the selection frame is displayed.
		
		<ul style="list-style-type: none"> With ▲▼ keys choose ACCELERATION to create acceleration characteristics or DECELERATION to create deceleration characteristics Confirm selection by pressing SET.

7.5.1 Acceleration characteristic, creation of episode 1

<p>Displayed alternately SPEED and 3000 (example):</p>  	No	section no. (max. 4)
	TIME	total acceleration time
	SPEED	final RPM
	ACC	characteristic's no. (10-19)
		adding a new section
		deleting last section
		editing sections
		exiting from characteristics wizard
		switching RPM/RCF
	entering graph view	

After entering the curve wizard, the symbol is highlighted. Pressing **SET** and selecting "NO" in response to the question "SAVE?" will return to the **PROG → CURVES** menu without making changes to the starting characteristics. To start editing the one-segment characteristics, select the icon with the ◀▶ keys and press the **SET** key.

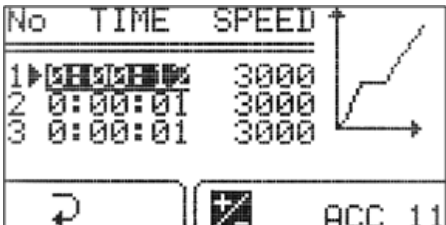
		editing value (flashing means editing the given value)
<ul style="list-style-type: none"> ▪ Press SET ▪ With ▲▼◀▶ choose time for section ▪ Press SET ▪ It is not possible to edit the maximum speed value. To do this, more sections have to be created, but the last section will always have the maximum set speed and cannot be changed. ▪ Select ↻ with ▼◀ buttons and press SET to finish editing characteristics. 		

7.5.2 Adding and editing sections - acceleration


To program next sections, select the ⊕ icon with the ◀▶ buttons and press **SET**. A new section (sections) will appear with a time of 1 second and a speed equal to the maximum speed.

To start editing a newly added section (sections), select the ✎ icon with the ◀▶ buttons and press **SET**, and follow the instructions given below.

After entering the profile section editing menu, the time value of the first section will be highlighted (see the picture below).

	The maximum speed value for the section cannot be higher than the maximum speed value for the characteristic (for the last section).
<ul style="list-style-type: none"> ▪ With ▲▼◀▶ highlight time or speed for desired section ▪ Press SET ▪ With ▲▼◀▶ choose value ▪ Press SET ▪ Repeat until setting all the sections ▪ To finish editing characteristic with ▲▼◀▶ choose ↻ and press SET. Finishing edition can be also done by pressing BACK button. 	


Saving created characteristic

	<ul style="list-style-type: none"> ▪ Select the ↻ icon with the ◀▶ buttons and press SET ▪ In the "Save?" window, use ▲▼ buttons to select YES to confirm saving the characteristic or NO, to exit without saving ▪ Press SET
---	--

7.5.3 Acceleration graph

An example of given parameters and a diagram:







No	TIME	9000
1	0:00:30	2000
2	0:00:52	6000
3	0:01:10	6000
4	0:02:00	9000



After the time value programming is completed, the **TIME + SPEED** segment of the user's startup characteristic can be displayed graphically. The set section of the characteristic curve is illustrated on the graph, which can be displayed by selecting the icon  with the ◀▶ keys and pressing the **SET** key.


7.5.4 Deceleration characteristic, creation of episode 1


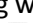

SPEED or **3000** displayed (example):

No	TIME	SPEED
1	0:00:15	0


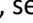

	adding a new section
	deleting last section
	editing sections
	exiting from characteristics menu
	switching RPM/RCF
	entering graph view


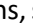

After entering the curve wizard, the symbol  is highlighted. Pressing **SET** and selecting "NO" in response to the question "SAVE?" will return to the **PROG → CURVES** menu without making changes to the starting characteristics. To start editing the one-segment characteristics, select the icon  with the ◀▶ keys and press the **SET** key.

	editing value (flashing means editing the given value)
<ul style="list-style-type: none"> ▪ Press SET ▪ With ▲▼◀▶ choose time for section ▪ Press SET ▪ To edit speed 	

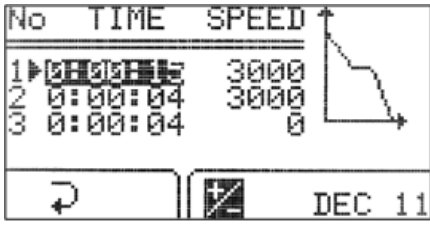












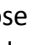
	<ul style="list-style-type: none"> It is not possible to edit the minimum speed value. To do this, more legends must be created, but the last leg will always be "0". Select  with   buttons and press SET to finish editing characteristics
--	--

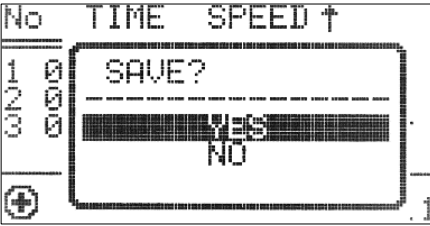




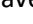
7.5.5 Adding and editing sections - deceleration

In order to program successive periods, select the icon  with the   keys and press the **SET** key. A new segment (or segments - after successive presses of SET) will appear with the time and speed equal to the minimum speed - "0".

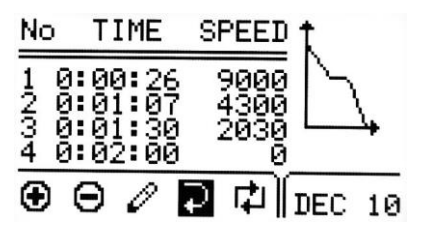

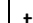

To start editing the newly added sections, select the icon  with the   buttons, press **SET** and make the settings as described below.

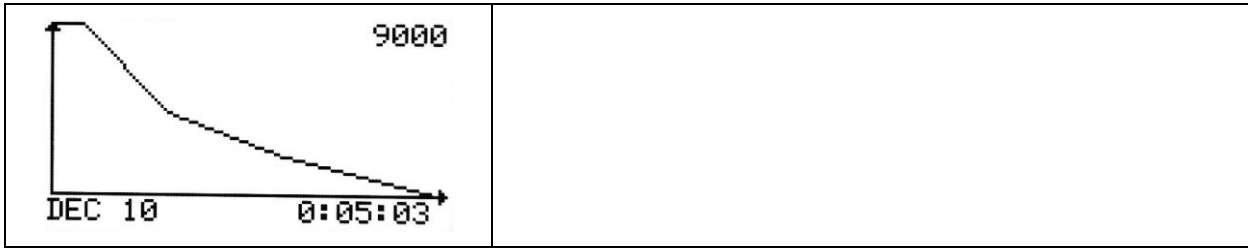
After entering the profile section editing menu, the time value of the first section will be highlighted (see the picture below).

	<p>The speed value of the last segment will always be "0".</p> <ul style="list-style-type: none"> With     highlight time or speed for desired section Press SET With     choose value Press SET Repeat until setting all the sections To finish editing characteristic with     choose  and press SET. Finishing edition can be also done by pressing BACK button
--	--

SAVING CREATED CHARACTERISTIC	
	<ul style="list-style-type: none"> Select the  icon with the   buttons and press SET In the "Save?" window, use   buttons to select YES to confirm saving the characteristic or NO, to exit without saving Press SET

7.5.6 Deceleration graph

An example of given parameters and a diagram:	
	<p>After programming the time value, you can graphically display the TIME + SPEED segment of the user's braking characteristic. The set section of the characteristic curve is illustrated on the graph, which can be displayed by selecting the icon  with the   keys and pressing the SET key.</p>



7.5.7 Deleting sections

In the characteristic's wizard:

	<ul style="list-style-type: none"> ▪ Select the icon with the buttons and press SET ▪ In the "Delete?" window, use buttons to select YES to confirm deleting the characteristic section or NO to cancel ▪ Press SET
--	---

7.6 Programs with user characteristics

Loading a modified program in the **CURVES** fold is signaled by the icon on the main screen:

	<p>Icon signals that program with user acceleration/deceleration characteristics are loaded.</p>
--	---

A change in any parameter entails the deactivation of the multi-section's curves mode.

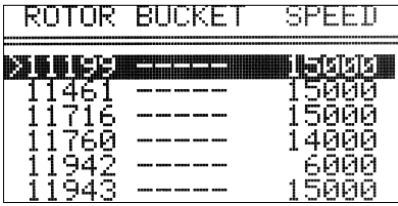
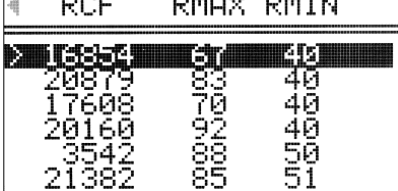
7.7 Rotor and bucket choosing

Simplified display mode


	<ul style="list-style-type: none"> ▪ Press and hold by 1 second. ▪ Choose rotor number (exemplary 11199/-----) with . ▪ Press SET. ▪ Execute points described follow (below Normal display mode description)
--	--

Normal display mode

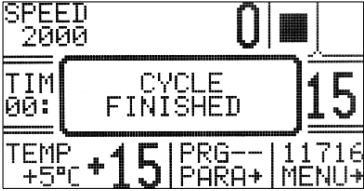


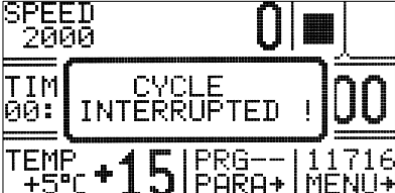
	<ul style="list-style-type: none"> ▪ Press SET— appears. ▪ Via mark rotor choosing field. ▪ Press SET (Rotors and buckets list will appear).
--	--

	<ul style="list-style-type: none"> Via ▲▼ keys mark demanded rotor number Confirm by pressing SET. If a bucket can be selected: <ul style="list-style-type: none"> With ▲▼ select demanded bucket number. Press SET. Press BACK to close edition mode.
	<ul style="list-style-type: none"> With ◀▶ keys one may switch between screens of rotors parameters
<ul style="list-style-type: none"> It is possible to set AUTOMATIC ROTOR IDENTIFICATION. The procedure is described in subsection <i>Other</i>. 	

7.8 SHORT mode

SHORT MODE – short work mode (centrifuging with pressed SHORT key)	
	<ul style="list-style-type: none"> The SHORT mode is activated by pressing and holding ►►(SHORT). In SHORT mode the centrifuge is working as long as the SHORT key is pressed or when set time is over. Centrifuging is stopped after releasing the SHORT key.

7.9 Finishing the centrifuging

<ul style="list-style-type: none"> When preselected time is reached, centrifugation will end automatically. 	
	
	<ul style="list-style-type: none"> Before lapsing preselected time, one may stop centrifugation. Pressing STOP for the first time will stop centrifuging with the characteristic set-in loaded program. ↓ symbol will be shown.
	<ul style="list-style-type: none"> Pressing STOP second time will stop centrifuging with the fastest characteristic. ↓ symbol will be shown.
	
<ul style="list-style-type: none"> The message about cancel of centrifuging can be delete with the STOP, SET, COVER, ▲▼◀▶ or BACK key. 	


7.10 Temporarily disabled functions

Functions written below can be temporarily disabled.





active	SPEED	RCF	TIME	TEMP	PROG —	— / — —	PARAM	MENU
THERMAL CHAMBER	●	●	●	○	●	●	●	●
STANDARD CENTRIFUGING	●	●	●	●	●	○	●	○

- available
- disabled


8. Temperature control

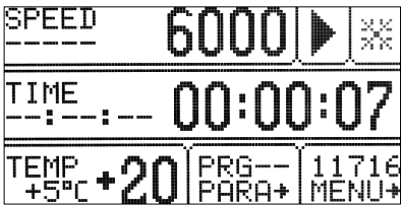

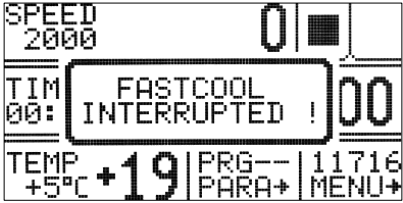
	MPW-260R only
---	----------------------

Centrifuge is equipped with ecological refrigerating system with temperature control. During centrifugation, there may appear differences in temperature on the display and temperature of the samples in the rotor. It depends on thermal conductivity of the rotor, and samples and centrifugation time, initial temperature of rotor and samples.


Exemplary change of TEMP setting:	
	<ul style="list-style-type: none"> ▪ Press SET (to enter edit mode)  – appears. ▪ Via ▲▼◀▶ keys mark TEMP field (highlighted). ▪ Press SET key. ▪ With ▲▼ choose demanded value (from -20°C to 40°C). ▪ Confirm settings by pressing SET. ▪ Press BACK.
	<p>Cooling is indicated by a symbol  (blinking).</p>

8.1 Initial cooling during centrifuging –FAST COOL


	<ul style="list-style-type: none"> ▪ The parameters allowable to change at FAST COOL mode: <ul style="list-style-type: none"> ▫ temperature (lower than current temperature shown by centrifuge) ▪ In order to centrifuging reduced temperature samples (e.g., storage in the external refrigerator) centrifuge chamber, rotor and centrifuge container must be pre-cooling to the predetermined temperature. It causes minimalization of temperature differences. ▪ Initial cooling may be activated by FAST COOL key (lid must be closed – rotor is spinning at FAST COOL mode) ▪ When FAST COOL mode is active, cooling system automatically set proper parameters to obtain demanded temperature the fastest way.
---	---

	<p>FAST COOL mode is marked by symbol  blinking in the right upper side of display.</p>
	<p>It is possible to exit FAST COOL mode at any time by pressing STOP key.</p> <p>Interruption of the function is signaled by a message.</p>


8.2 Initial cooling without centrifuging – THERMAL CHAMBER

	<p style="text-align: right;">PARAM → THERMAL CHAMBER</p> <ul style="list-style-type: none"> There is possible to run centrifuge in THERMAL CHAMBER mode – cooling (rotor is at standstill). How to enable THERMAL CHAMBER is described in “Thermal chamber” chapter.
---	--

8.3 Cooling in “START DELAY – OF TEMPERATURE” mode

	<p style="text-align: right;">PARAM → START DELAY – OF TEMPERATURE</p> <ul style="list-style-type: none"> Centrifuging process will start, when preselected temperature is reached. How to enable run START DELAY – OF TEMPERATURE function is described in Start delay – of temperature chapter.
---	--

8.4 Cooling in „SHORT” mode



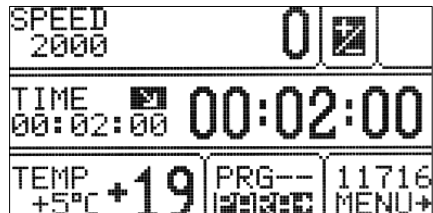


	<ul style="list-style-type: none"> Cooling feature is available in SHORT mode. How to enable run centrifugation in SHORT mode is described in SHORT mode.
---	--

8.5 Cooling notes

Centrifuge with cooling – MPW-260R is equipped with an efficient cooling system. It allows obtaining selected temperatures in the chamber even at maximum spin speed or fast obtaining desired temperatures (e.g., +4°C). Note that time and possibility of obtaining a set temperature is dependent on multiple factors, including: the power of the cooling system, the shape of the rotor, the rotor speed, ambient temperature, etc. The temperature on the display is appropriate for the place of the temperature sensor in the chamber, accuracy is ±1°C. The temperature of the sample may be different.


9. Parameters of centrifugation

This chapter contains exemplary screens of MPW-260R centrifuge (screens for MPW-260 – without cooling – do not include temperature field).

Simplified display	
	<ul style="list-style-type: none"> Press and hold  by 1 second. Choose PARAM. with ▲▼ Press SET.
Normal display	
	<ul style="list-style-type: none"> Press SET. With ▲▼◀▶ keys select PARAM. Press SET.
<ul style="list-style-type: none"> It is possible to switch between two different screens via ◀▶ keys in PARAMETERS field. 	
	

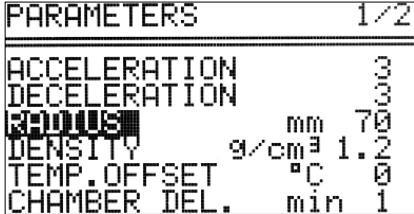
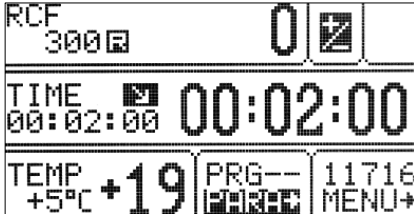
ACCELERATION	chosen acc. characteristic (0-the fastest, 9-the slowest)
DECELERATION	chosen dec. characteristic (0-the fastest, 9-the slowest)
RADIUS [mm]	current rotor radius [mm]
DENSITY (g/cm³)	sample density [g/cm ³]
TEMP. OFFSET (°C)	value of temperature correction
CHAMBER DEL. (min)	delay between set thermal chamber mode and start it
THERMAL CHAMBER	cooling of the chamber without centrifuging
AUTOM. LID OPENING	opening cover after centrifuging automatically
START DELAY	starting delayed (after pressing START)

9.1 Acceleration/deceleration – changing characteristics

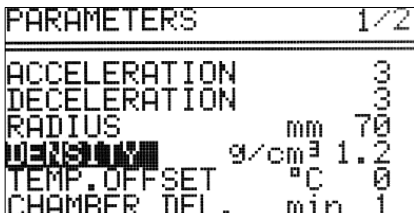
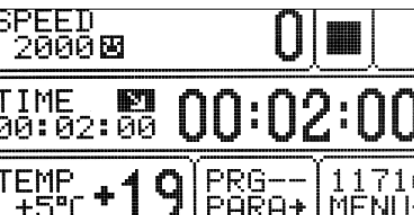
	<ul style="list-style-type: none"> With ▲▼ keys select ACCELERATION or DECELERATION. Press SET. With ▲▼ keys select demanded number of characteristics. Press SET.
---	--

	<p>ACCELERATION –10 (0 ÷ 9), linear accelerating characteristics assigned to every rotor. 0-the fastest acceleration, 9-the slowest acceleration.</p> <p>DECELERATION – 10 (0 ÷ 9), linear decelerating characteristics assigned to every rotor. 0-the fastest deceleration, 9-the slowest deceleration.</p>
--	--

9.2 Radius

 <p>PARAMETERS 1/2</p> <p>ACCELERATION 3</p> <p>DECELERATION 3</p> <p>RADIUS mm 70</p> <p>DENSITY g/cm³ 1.2</p> <p>TEMP.OFFSET °C 0</p> <p>CHAMBER DEL. min 1</p>	<p>RADIUS [mm] - control of the radius of the rotor within the range from R_{min} to R_{max}. Available values depend on chosen rotor. Radius corrections serve for more precise control RCF, exemplary when user need to know real RCF in half length of test tube.</p> <ul style="list-style-type: none"> To change the rotor radius, select RADIUS [mm] with ▲▼ keys. Press SET. Set demanded value by pressing ▲▼. Press SET.
 <p>RCF 300 R 0 R</p> <p>TIME 00:02:00</p> <p>TEMP +5°C +19 PRG-- 11716 MENU+</p>	<p>When radius correction is activated, R symbol is visible on the screen.</p> <p>Reducing of the rotor radius resulting change of displayed RCF value.</p>

9.3 Sample density

 <p>PARAMETERS 1/2</p> <p>ACCELERATION 3</p> <p>DECELERATION 3</p> <p>RADIUS mm 70</p> <p>DENSITY g/cm³ 1.2</p> <p>TEMP.OFFSET °C 0</p> <p>CHAMBER DEL. min 1</p>	<p>DENSITY (g/cm³) – default density is set to 1,2 g/cm³</p> <p>To change the density (possible values 1,2÷9,9 g/cm³):</p> <ul style="list-style-type: none"> Via ▲▼ keys select DENSITY (g/cm³) Press SET. Set demanded value by pressing ▲▼. Press SET.
 <p>SPEED 2000 D 0 D</p> <p>TIME 00:02:00</p> <p>TEMP +5°C +19 PRG-- 11716 MENU+</p>	<p>When density is changed, D symbol is visible on the screen.</p> <p>Changing of DENSITY value is obligatory when density of sample placed into rotor is higher than 1.2 g/cm³. Change of DENSITY value led to decreasing maximum value of accessible speed.</p>

9.4 Temperature offset



MPW-260R only

```

PARAMETERS 1/2
-----
ACCELERATION 3
DECELERATION 300
RADIUS mm 70
DENSITY g/cm³ 1.2
TEMP.OFFSET °C 0
CHAMBER DEL. min 1
    
```

Temperature offsets serve for more precise control of real sample temperature. It can be helpful in case high/low initial temperature samples or high-volume samples.

- With ▲▼ keys select **TEMP. OFFSET**.
- Press **SET**.
- Use the ▲▼ keys to select the difference between the temperature that the cooling system will aim for and set temperature. Confirm selection by pressing **SET**.

Attention!

The use of the offset cannot extend the temperature range achieved by the centrifuge.

Function description

At a set temperature of 20°C and the set offset value equal to -5°C, cooling system will actually strive to reach 15°C. With a setpoint temperature of 20°C and a set offset value of 5°C the system will actually try to reach 25°C.



The temperature displayed on the main screen is corrected for offset value.

Offset can be selected range from -20°C to 20°C.

```

SPEED 176
2000
-----
TIME 00 01 59
00:02:00
-----
TEMP +5°C +21 PRG-- 11716
+5°C +21 PARA+ MENU+
    
```

Activation of the function is signaled on the main screen

with  or  depending on the offset value sign.

9.5 Thermal chamber



MPW-260R only

Cooling without centrifuging.

THERMAL CHAMBER

```

PARAMETERS 2/2
-----
[ ] THERMAL CHAMBER
[ ] AUTOM. LID OPENING
[ ] START DELAY
    
```

- With ▲▼◀▶ keys select **THERMAL CHAMBER**.
- Press **SET** (to turn on/off).
- With ▲▼ keys select temperature value.
- Set demanded value (0°C÷40°C) by pressing ▲▼.
- Confirm selection by pressing **SET**.

Attention, in the centrifuge without heating, do not set the thermal chamber to a value higher than currently indicated by the centrifuge.


	<p>When THERMAL CHAMBER function is activated, T symbol is visible on the screen.</p> <p>Changing temperature from the main screen is not possible.</p> <p>Opening cover terminates THERMAL CHAMBER function (closing cover back turns it on).</p>
	<p>Thermal chamber is activated with delay.</p> <ul style="list-style-type: none"> Set time of delaying by select CHAMBER DEL. Press SET. With ▲▼ keys select demanded value (1-5 min). Press SET.
<ul style="list-style-type: none"> The function is activated automatically after confirmation and with the lid closed. The function is interrupted when the lid is opened, and the function resumes when the lid is closed again. If the THERMAL CHAMBER function is enabled during the centrifugation cycle, at the end of this cycle, the THERMAL CHAMBER function is activated until the lid is opened. Unlike other parameters, the THERMAL CHAMBER function can be turned on only when the centrifuge is stopped. 	

9.6 Automatic lid opening

<p>Automatic lid opening</p>	<p style="text-align: center;">AUTOMATIC LID OPENING</p>
	<ul style="list-style-type: none"> When centrifuge process is finished, cover will be opened automatically for set option AUTOM. LID OPENING. When centrifuging is terminated by pressing STOP, opening cover is possible by pressing COVER.
	<p> symbol means that OPEN LID AFTER RUN is active.</p>


9.7 Start delay - of time

	<p>Start centrifuging since preselected delay is reached.</p>	<p style="text-align: center;">START DELAY / OF TIME</p>
		<ul style="list-style-type: none"> With ▲▼ keys select START DELAY. Press SET. Start delay can be set from 0 : 0 0 : 0 1 to 9 : 5 9 : 5 9. With ▲▼ keys select OF TIME. Press SET and ▶ and then SET. With ▲▼ keys set demanded value. Confirm by pressing SET. Press BACK to escape edit mode.

	<p>When START DELAY function is activated,  symbol is visible on the screen.</p>
<ul style="list-style-type: none"> START DELAY / OF TIME function can be stopped at any moment by pressing STOP. 	
<ul style="list-style-type: none"> START DELAY / OF TIME function cannot be run when START DELAY / OF TEMP. is activated. 	

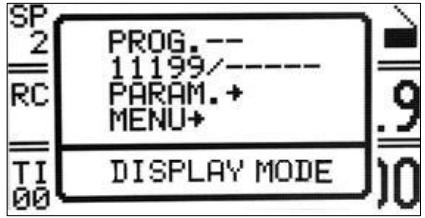



9.8 Start delay – of temperature

	<p>MPW-260R only</p>
--	----------------------


	<p>Start centrifuging time counting since preselected temperature is reached.</p>	<p>START DELAY / OF TEMP.</p>
		<ul style="list-style-type: none"> With ▲▼◀▶ keys mark START DELAY. Press SET. With ▲▼◀▶ keys mark OF TEMP. Press SET. Press ▶, press SET. With ▲▼ keys set demanded value of temperature. Press SET. Exit edit mode by press BACK.
		<p>When START DELAY – OF TEMPERATURE is turned on,  symbol is visible on the screen.</p>
<ul style="list-style-type: none"> When the function is active, the speed can be reduced to the optimum values for the FAST COOL function, when the set speed is lower than the optimum value, the rotor rotates at the set by user speed. 		
<ul style="list-style-type: none"> The delay starts - from the temperature can be interrupted at any time by pressing the STOP key. 		
<ul style="list-style-type: none"> START DELAY / OF TEMP. function cannot be run when START DELAY / OF TIME is activated. 		

10. Menu

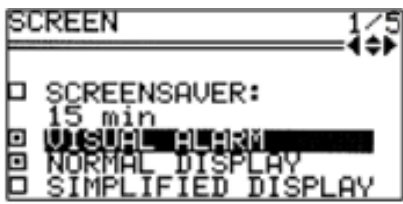
This chapter contains exemplary screens of MPW-260R centrifuge (screens for MPW-260 – without cooling – do not include temperature field).

Simplified display	
	<ul style="list-style-type: none"> Press and hold  by 1 second. Choose MENU with ▲▼ Press SET. <p>Execute points described follow (below Normal display mode description)</p>
Normal display	
	<ul style="list-style-type: none"> Press SET. With ▲▼◀▶ keys select MENU. Press SET.
	<ul style="list-style-type: none"> To navigate in MENU use ▲▼◀▶ keys. To enter menu press SET.
CONFIGURATION	centrifuge configuration
PASSWORD	password protection
LAST 10-CYCLES	10 last centrifugation cycles history
CYCLES	total working time of centrifuge, total number of working cycles
ROTOR RUNTIME	counting time of work and cycles amount for each rotor
CONTACT US	manufacturer's details
DIAGNOSTICS	error codes (service field)
FACTORY SETTINGS	restore factory settings

10.1 Screen saver

Setting time of screen saver	MENU/ CONFIGURATION / SCREEN
	<ul style="list-style-type: none"> With ▲▼◀▶ keys select SCREENSAVER. Press SET and then ▼ and SET. With ▲▼ keys select demanded value from 1 to 60 minutes. Mark selection by pressing SET. Leave the menu by pressing BACK.

10.2 Visual alarm

Visual alarm	MENU/CONFIGURATION/ SCREEN
	<ul style="list-style-type: none"> Via ▲▼ keys choose VISUAL ALARM Mark it by pressing SET. <p>VISUAL ALARM cause blinking screen after ending of centrifuging or after message occurring.</p>

10.3 Types of main screen

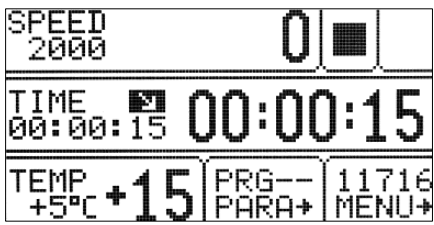
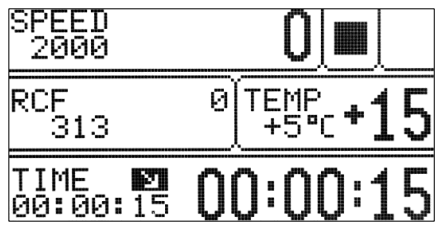
To ensure optimal adaptation to the user's preferences, work is possible in two basic screen modes.

NORMAL DISPLAY - contains an expanded number of parameters visible on the display.

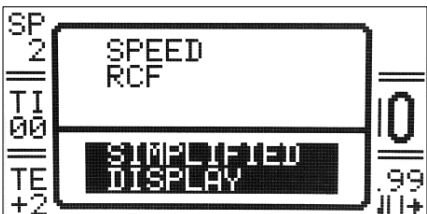
SIMPLIFIED DISPLAY - contains only the most important parameters visible on the display.

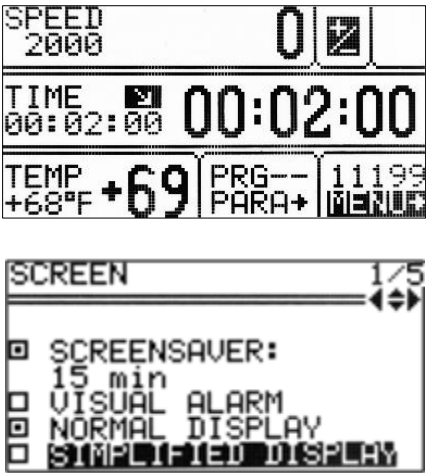

For each of the above modes, you can choose priority RPM display or RCF.

By default, the **SIMPLIFIED DISPLAY** is set

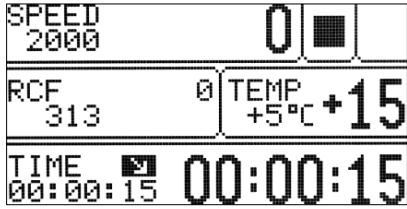
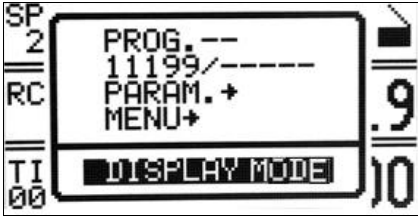
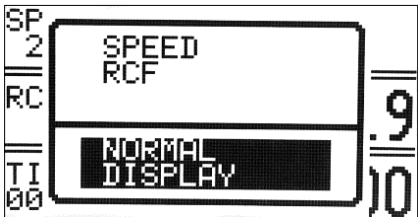
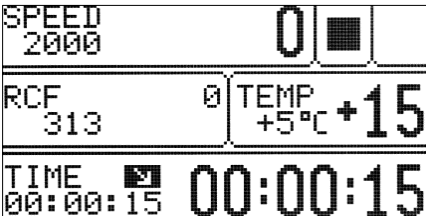
Types of main screen	
NORMAL DISPLAY	SIMPLIFIED DISPLAY
	
Switch between the SPEED (RPM) and RCF display priority modes	
<ul style="list-style-type: none"> In the NORMAL DISPLAY mode, selecting the SPEED or RCF display mode is obtained by pressing and holding BACK for 1 sec. then use the ▲▼ buttons to select the desired mode (SPEED or RCF) and press SET. 	<ul style="list-style-type: none"> In the SIMPLIFIED DISPLAY mode, the selection of the SPEED or RCF display mode is obtained by pressing and holding the BACK key for 1 second. then use ▲▼ keys to select DISPLAY MODE, press SET, and then use ▲▼ keys to select the desired mode (SPEED or RCF) and press SET.

10.3.1 Switching the normal display to simplified display

Method I	
	<ul style="list-style-type: none"> Press the BACK button for 1 sec. to return to the basic display (a short menu is displayed on the screen), then: Via ▲▼ keys select SIMPLIFIED DISPLAY. Press SET.

Method II	
	<ul style="list-style-type: none"> ▪ Press SET –  appears. ▪ Via ▲▼◀▶ keys select MENU. ▪ Press SET. ▪ Via ▲▼ keys select CONFIGURATION tab. ▪ Press SET. ▪ Via ◀▶ keys select SCREEN tab. ▪ Via ▲▼ keys select SIMPLIFIED DISPLAY. ▪ Press SET. ▪ Leave menu via BACK key.

10.3.2 Switching the simplified screen to normal display

Method I	
	<ul style="list-style-type: none"> ▪ Press the BACK button for 1 sec.
	<ul style="list-style-type: none"> ▪ Via ▲▼ keys select DISPLAY MODE (highlighted). ▪ Press SET.
	<ul style="list-style-type: none"> ▪ Then choose NORMAL DISPLAY via ▲▼ keys. ▪ Press SET.
Method II	
	<ul style="list-style-type: none"> ▪ Press the BACK button for 1 sec.

	<ul style="list-style-type: none"> Via ▲▼ keys select MENU (highlighted). Press SET.
	<ul style="list-style-type: none"> Via ▲▼ keys select CONFIGURATION tab. Press SET. Via ◀▶ keys select SCREEN tab. Via ▲▼ keys select NORMAL DISPLAY. Press SET. Leave menu via BACK key.

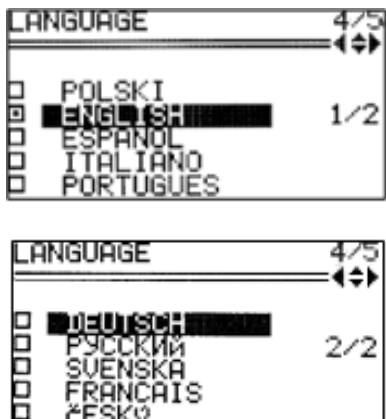
10.4 Rotating runtime

Way of time counting	MENU/CONFIGURATION/ ROTATING RUNTIME
	<ul style="list-style-type: none"> Via ▲▼ choose demanded option. Mark it by pressing SET.
Counting from:	
From pressing start →	COUNTING SINCE ROTOR IS IDENTIFIED
From reaching speed →	COUNTING FROM ASSIGNED SPEED
Presenting mode:	
Descending →	COUNTING DOWN
Ascending →	COUNTING UP

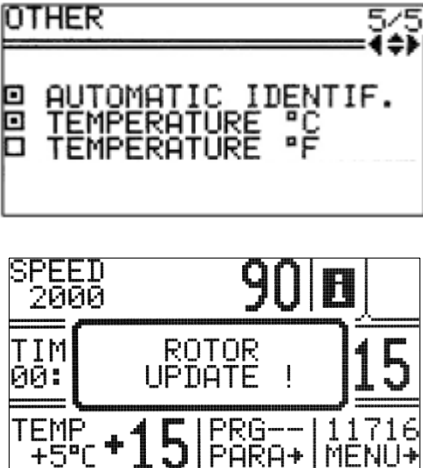
10.5 Buzzer

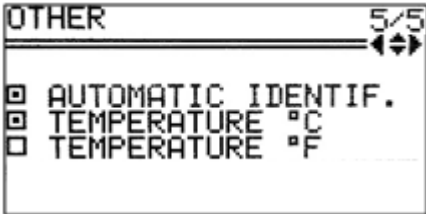
Switching ON/OFF short audible signals accompanying every pressing of any key. Switching ON/OFF signals after centrifuging.	MENU/ CONFIGURATION /BUZZER
	<ul style="list-style-type: none"> With ▲▼ keys select demanded option. Mark selection by pressing SET. <p>A continuous alarm means the emission of short beeps after the end of the spin, until the message about the end of the work cycle is deleted.</p>
<ul style="list-style-type: none"> Warning signals are always switched on. 	

10.6 Language

Changing menu language	MENU / CONFIGURATION / LANGUAGE
	<ul style="list-style-type: none"> ▪ Via ▲▼ keys choose demanded menu language ▪ Mark it by pressing SET.

10.7 Other

Rotor automatic identification	MENU / CONFIGURATION / OTHER
	<p>Thanks to the automatic rotor identification, the centrifuge automatically identifies the rotor in the chamber. Rotor identification is indicated by the message.</p> <p>When the function is deactivated, it is necessary to manually select the desired rotor as described in “Rotor and bucket choosing”.</p> <p>The AUTOMATIC IDENTIF. is turned on by default.</p> <p>To enable/unable the function:</p> <ul style="list-style-type: none"> ▪ Via ▲▼ keys choose <ul style="list-style-type: none"> <input type="checkbox"/> AUTOMATIC IDENTIF. ▪ Press SET (<input type="checkbox"/> change to <input checked="" type="checkbox"/> or conversely). <p>Autoidentification is not active for work in the loaded program mode.</p>

Choice of temperature unit (only MPW-260R)	MENU / CONFIGURATION / OTHER
	<p>The TEMPERATURE in °C is turned on by default.</p> <p>To change the temperature unit:</p> <ul style="list-style-type: none"> ▪ Via ▲▼ keys select unit ▪ Confirm by pressing SET.

--	--

10.8 Password

Setting up password	MENU / PASSWORD
---------------------	-----------------

To prevent from an unauthorized use, a **PASSWORD** can be set.

Note: No PASSWORD is set by default.

The PASSWORD can be set as follows when the rotor is at a standstill.

	<ul style="list-style-type: none"> Press SET. Icon starts blinking. With ◀▶ keys set the valid place of the PASSWORD. With ▲▼ keys set correct value. Repeat above steps for all places. Press SET.

	<ul style="list-style-type: none"> As a confirmation repeat instructions described above.
--	--

When the **PASSWORD** is set, the Key sign is displayed in the **CODE** zone. It is also displayed in the main menu (lower right corner of the screen).

--	--

- From then on, access to the **MENU** is possible after entering the password.
- In case of incorrect password, it will show message: **ACCESS DENIED!**

- To delete the **PASSWORD**, "0000" must be set (after previously entering current password).

If the **PASSWORD** is forgotten, the emergency code “**7654**” should be used to clear password and remove all locks.

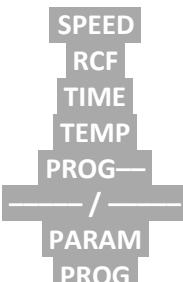
Setting up locks



- With ▲▼ keys choose a lock.
- Mark a lock by pressing **SET**.
- Repeat above steps for desired locks.
- Leave menu with **BACK** key.



- In the LOCK menu, press ◀ and then SET.

	disabled*	description
SAVE PROGRAM	SAVE button	no programs can be saved
DELETE PROGRAM	DELETE button	no programs can be deleted saving programs on position where one was already stored is disabled
CHANGE PARAMETERS	fields: 	parameters cannot be modified
LOAD PROGRAM	LOAD button	no programs can be called up
START KEY	START key	centrifugation cannot be started

* Executing disabled procedures is only possible after entering the correct password.

10.9 Last 10 cycles

Information concerning parameters of last 10 centrifuging cycles.

MENU / LAST 10 CYCLES



- Number of cycles can be changed by ◀▶ keys.
- The list can be scrolled using ▲▼ keys.
- To exit press **SET/BACK** key

10.10 Work time

<p>Total working time of centrifuge, and quantity of working cycles.</p>	<p>MENU / WORK TIME</p>
<pre> WORK TIME ----- TOTAL RUN TIME: 0h 13m 14s CYCLES: 31 </pre>	<p>In the WORK TIME menu, the following statistics are displayed:</p> <ul style="list-style-type: none"> ▪ total working (centrifugation) time ▪ working cycles counter

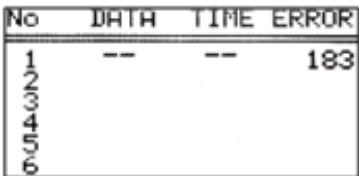
10.11 Rotor runtime

<p>Information about the time of centrifuging and of the quantity of the working cycles of each rotor. The table also contains icons warning of the duty of execution of validation.</p>	<p>MENU / ROTOR RUNTIME</p>
<pre> ▶ ROTOR CYCLES NOM.C. ----- ✓ 11199 0 15000 ✓ 11461 0 15000 ✓ 11716 15 15000 ✓ 11760 0 15000 ✓ 11942 11 15000 ✓ 11943 0 15000 </pre>	<p>CYCLES – the number of centrifuging the rotor has performed, NOM.C. – permissible number of centrifuging for the rotor.</p> <ul style="list-style-type: none"> ▪ The list can be scrolled using ▲▼ keys. ▪ To exit press BACK key. <p>Symbols:</p> <ul style="list-style-type: none"> ✓ – more than 100 cycles left ! – less than 100 cycles left ■ – worn rotor <p>It is not allowed to use rotors marked as worn.</p>

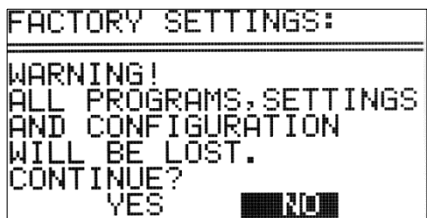
10.12 Contact us

<p>Information about the type of the centrifuge, firmware version, and contact details.</p>	<p>MENU / CONTACT US</p>
<pre> CONTACT US ----- MPW MED. INSTRUMENTS 04-347 WARSAW 46 BOREMLOWSKA Street ----- WWW.MPW.PL MPW@MPW.PL </pre>	<ul style="list-style-type: none"> ▪ The list can be scrolled using ▲▼ keys. ▪ To exit press BACK key.

10.13 Diagnostics


Information about errors arisen in working of the centrifuge (for service).	MENU / DIAGNOSTICS
	Intended for service purposes!

10.14 Factory settings

Restoring factory settings.	MENU/ FACTORY SETTINGS
All settings of user programs will be deleted.	
	<ul style="list-style-type: none"> ▪ Via ◀▶ keys choose YES or NO. ▪ Confirm by pressing SET.

11. Maintenance

11.1 Cleaning of the centrifuge

	<ul style="list-style-type: none"> ▪ Pull the mains plug before cleaning. ▪ Before any cleaning or decontamination process other than that is recommended by the manufacturer, the user has to ask the manufacturer if the planned process does not damage the device ▪ For cleaning, water with soap or other water-soluble mild detergent shall be used. ▪ One should avoid corrosive and aggressive substances. It is prohibited to use alkaline solutions, inflammable solvents or agents containing abrasive particles. ▪ Do not lubricate the centrifuge motor shaft. ▪ The unused centrifuge should have cover opened. Once a week ▪ Using wiping cloth, remove condensate or residues of the products from the rotor chamber. Once a month ▪ Check the rotor fixing screw thread. In case of damage, replaced it. ▪ Check the centrifuging chamber whether it is damaged. In case of damage, it cannot be longer put into operation. Notify authorized service workshop.
---	--

11.2 Maintenance of centrifuge elements



The rotor pins shall be always lubricated with petroleum jelly. In this way, the uniform deflection of the buckets and quiet centrifuge operation is ensured.

Cleaning of the accessories



- In order to ensure safe operation, one shall carry out in **regular** way periodical maintenance of the accessories.
- Rotors, buckets and round carriers have to withstand high stresses originating from the centrifugal force. Chemical reactions as well as corrosion (combination of variable pressure and chemical reactions) can cause destruction of metals. Hard to observe surface cracks increase gradually and weaken material without visible symptoms.
- In case of observation of surface damage, crevice or other change, as well as the corrosion, the given part (rotor, bucket, etc.) shall be immediately replaced.
- The rotor, including the fixing screw, buckets and round carriers must be regularly cleaned to prevent corrosion.
- Cleaning of the accessories shall be carried out outside of the centrifuge **once every week** or still better after each use. For cleaning them one should use neutral agent of pH value 6÷8. It is forbidden to use alkaline agent of **pH > 8**. Then, those parts shall be dried using soft fabric or in the chamber drier at ca. 50°C.
- Angle rotor should be placed on a fabric with holes facing down, for effective drying.
- Do not use bleach on plastic parts of the rotor.
- In this way, the useful service life of the device is substantially increased and susceptibility to corrosion is diminished. Accurate maintenance increases the service life as well and protects against premature rotor failures.
- Do not use bleach on plastic parts of the rotor.
- According to laboratory standards, minimize the immersion time in each solution.
- Especially prone to the corrosion are parts made of aluminum.
- Corrosion and damages resulting from insufficient maintenance could not be subject of claims lodged against the manufacturer.
- The unused rotor should have the lid removed.

HS accessories maintenance.



- Check the general condition of seals.
- Make sure that rubber O-rings are lightly coated with silicone grease. Use high vacuum grease, e.g., type „C“ by LUBRINA.
- In order to maintain hermetic sealing, it is recommended to replace the sealing rings after each autoclaving.
- Store hermetically sealed rotors and buckets with the lids removed.

11.3 Sterilization

Plastics - legend to abbreviations

PS	polystyrene	ECTFE	ethylene/chlorotrifluoroethylene
SAN	styrene-acrylonitrile	ETFE	ethylene/tetrafluoroethylene
PMMA	polymethyl methacrylate	PTFE	polytetrafluoroethylene
PC	polycarbonate	FEP	tetrafluoroethylene/perfluoro propylene
PVC	polyvinyl chloride	PFA	tetrafluoroethylene/perfluoroalkylvinylether

POM	acetal polyoxymethylene	FKM	fluorocarbon rubber
PE-LD	low density polyethylene	EPDM	ethylene propylene diene
PE-HD	high density polyethylene	NR	natural rubber
PP	polypropylene	SI	silicon rubber
PMP	polymethyl pentene		

One can use all standard disinfectants. Centrifuges and devices are made of different materials, one should consider their variety.

	radiation β radiation γ 25 kGy	C ₂ H ₄ O (ethylene oxide)	formalin, ethanol
PS	●	○	●
SAN	○	●	●
PMMA	●	○	●
PC	●	●	●
PVC	○	●	●
POM	●	●	●
PE-LD	●	●	●
PE-HD	●	●	●
PP	●	●	●
PMP	●	●	●
ECTFE, ETFE	○	●	●
PTFE	○	●	●
FEP, PFA	○	●	●
FKM	○	●	●
EPDM	○	●	●
NR	○	●	●
SI	○	●	●

- may be used
- cannot be used

In the centrifuge, disinfectants and cleaning agents generally used in medical care should be used (e.g., *Aerodesina-2000*, *Lysoformin 3000*, *Melseptol*, *Melsept SF*, *Sanepidex*, *Cutasept F*).

11.3.1 Autoclaving

- Rotors, buckets and round carriers can be sterilized in autoclave with temperature 121°C during 20 min (215 kPa), unless otherwise specified in the OPTIONAL ACCESSORY.
- During sterilization (autoclaved) by means of steam one should consider temperature resistance of individual materials.
- Deformation of the accessories (carriers or lids made of plastic) may occur during autoclaving.
- Do not autoclave disposable materials (e.g., tubes, cyto-container).
- The life of the accessory depends on the frequency of autoclaving and use.
- Autoclaving reduces lifespan of plastic components. They should be replaced if any signs of damage are visible, including a change in color or shape or when leakage etc.
- Pressure in closed containers can cause plastic deformation or explosion.
- Prior to autoclaving the rotors and accessories, thoroughly wash and rinse with distilled water.
- Never exceed the permissible autoclaving temperature and time.
- If you want to keep the hermetic seals, replace the sealing rings after each autoclave.

Chemical resistance of plastics

	autoclaving 121 °C, 20 min		autoclaving 121 °C, 20 min
PS	○	PMP	●
SAN	○	ECTFE, ETFE	●
PMMA	○	PTFE	●
PC	●	FEP, PFA	●
PVC	○ ¹⁾	FKM	●
POM	●	EPDM	●
PE-LD	○	NR	○
PE-HD	○	SI	●
PP	●		

● may be used

○ cannot be used

1) Except PVC hoses which are resistant to the steam sterilization in the temperature 121°C.

11.4 Chemical resistance

Chemical resistance of plastics

	aldehydes	cyclic alcohols	esters	ether	ketones	strong or concentrated acids	weak or diluted acids	oxidizing substances	cyclic hydrocarbons	aromatic hydrocarbons	haloid hydrocarbons	alkalis
PS	○	●	○	○	○	○/●	○/●	○	○	○	○	●
SAN	○	●	○	○	○	○	○/●	○	○	○	○	●
PMMA	○/●	●	○	○	○	○	○/●	○	○/●	○	○	○
PC	○/●	●	○	○	○	○	○/●	○	○/●	○	○	○
PVC	○	●	○	○	○	●	●	○	●	○	○	●
POM	○/●	●	○	●	●	○	○	○	●	●	●	●
PE-LD		●	●	●	○/●	●	●	○	●	●	●	●
PE-HD	●	●	○/●	○/●	○/●	●	●	○	●	○/●	○/●	●
PP	●	●	○/●	○/●	○/●	●	●	○	●	○/●	○/●	●
PMP	○/●	●	○/●		○/●	●	●	○	○/●	○	○	●
ECTFE ETFE	●	●	●	●	○	●	●	●	●	●	●	●
PTFE FEP PFA	●	●	●	●	●	●	●	●	●	●	●	●
FKM	●	○	○	○	○	○	●	○/●	○/●	○/●	○/●	○/●
EPDM	●	●	○/●	○	○/●	●	●	○/●	○	○	○	●
NR	○/●	●	○/●	○	○	○	○/●	○	○	○	○	●
SI	○/●	●	○/●	○	○	○	○/●	○	○	○	○	○/●

● very good

Permanent action of the substance does not cause damage through 30 days. The material is able to be resistant through years

○/● good to limited



Continuous action of the substance causes insignificant and partly reversible damage through the period of 7-30 days (e.g., puffing up, softening, reduced mechanical durability, discoloring).

○ limited

The material should not have the continuous contact with the substance. The immediate occurrence of damage is possible (e.g., the loss of mechanical durability, deformation, discoloring, bursting, dissolving).

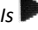

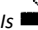

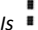
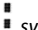
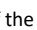
Rubber inserts shall be exactly cleaned or possibly replaced. Centrifuges and accessories are made of different materials.

Do not use bleach on plastic parts of the rotor.

	DANGER! MPW accessories are not biotight. For centrifuging infectious materials, it is necessary to use hermetically closed tubes meeting demands of biotightness, in order to prevent germs migration into the centrifuge and beyond it.
	User is responsible for proper disinfections of the centrifuge if some dangerous material was spilled inside or outside of the centrifuge. During the above-mentioned works, one must wear safety gloves.

12. Troubleshooting

Majority of faults could be removed by switching the centrifuge OFF and then ON. After switching the centrifuge ON, there shall be displayed parameters of the recently implemented program and sound signals comprising four successive tones shall be generated. In case of short-duration power failure the centrifuge terminates the cycle and displays PROGRAM ERROR code.

problem	question	remedy
Centrifuge does not start	<i>Is supply cable plugged into mains?</i>	<i>Plugs supply cable correctly.</i>
	<i>Is master switch ON?</i>	<i>Switch ON power supply.</i>
Motor error is displayed		Call service.
Centrifuge does not start (indications are proof for cycle in progress and motor does not start)	<i>Is  symbol displayed?</i>	Wait till rotor stops and the  symbol goes off.
	<i>Is  symbol displayed?</i>	Close cover.  symbol must switch off.
	<i>Is  symbol blinking?</i>	Centrifugation cycle in progress, press STOP key or wait till cycle ends.
Centrifuge does not accelerate (unbalance error)	<i>Unequal rotor load.</i>	Centrifuge load shall be balanced.
	<i>Inclined centrifuge.</i>	Centrifuge shall be levelled.
	<i>Faulty drive (mechanical damage).</i>	Call service.
	<i>Was centrifuge displaced during operation?</i>	Switch ON the centrifuge again after opening and closing the cover.
(motor error)	<i>After stopping error rotor message is displayed</i>	Check if rotor number in started program is consistent with the number of the rotor installed in the centrifuge. Check rotor status (if there are coding magnets inserted)
	<i>Centrifuge does not recognize the rotor and does not stop.</i>	Switch the centrifuge OFF, then ON and check correctness of loaded program
It is not possible to open the cover	<i> symbol on the display is blinking, after pressing COVER key single tone is audible</i>	Rotor is still rotating. Wait for stopping of the rotor and displaying of the  symbol.
	<i>The sensor is connected correctly, and the error is still applying.</i>	Call service.
Mains failure during run	<i>The message will be displayed on the display about the decay of tension.</i>	Wait for stopping of the rotor, clear the error by pressing the SET key.
Temperature sensor error	<i>The overheating message will be displayed.</i>	Switch the centrifuge OFF, then ON.
		Call service.


Error of the exceeding the temperature (50°C) in the chamber	<i>The overheating message will be displayed.</i>	Open the cover. Wait for the centrifuge to cool down.
--	---	---

12.1 Messages

Screen messages that may occur during operation.	
MESSAGE	EXPLANATION
"SPEED OF ROTOR" "IDENTIFICATION <> 90 RPM"	Please try start centrifuging again, if error still occur, contact manufacturer's authorized service.
"IMBALANCE FAST STOP !" "PLEASE REMOVE CAUSE" "THEN RESTART"	Rotor is not balanced correctly, please balance rotor.
"NO ROTOR OR IDENTIFICATION" "SENSOR DAMAGED !"	Make sure, is rotor mounted in the centrifuge chamber. If it is right contact manufacturer's authorized service.
"INCORRECT ROTOR NUMBER !"	Change rotor number in centrifuge settings or use autoidentification.
"WRONG DIRECTION OF ROTATION" "OR UNKNOWN ROTOR !"	Make sure if correct rotor for centrifuge is mounted. List of accessories is described in chapter 15.
"PLEASE CLOSE THE LID" "HAND !"	Necessity of manually closing the lid.
"ROTOR STOPPING !" "Please wait..."	Initializing after mains failure with rotating rotor, wait until rotor stop.

Emergency messages
In case of emergency messages (centrifuge is not working properly) contact the manufacturer's authorized service center.
MESSAGE
"OVERHEATING MOTOR !" "INVERTER ERROR !"
"INVERTER SERIAL BUS ERROR !"
"TEMPERATURE SENSOR ERROR"
"PRESSURE CONTROL FAILURE!"
"OPENING COVER in RUN!"
"SPEED METER ERROR"
"I2C BUS ERROR"
"OVERHEATING CENTRIFUGE !"
"ROTOR OVERSPEED !"
"COVER LOCK MALFUNCTION !"

12.2 Emergency cover release


	<p>EMERGENCY COVER RELEASE</p> <p>Attention! <i>The cover may be opened in emergency only when the rotor is at rest. Before emergency opening the cover, switch off the mains power switch and disconnect the power cord. Wait 10 min and/or looking through the sight glass, make sure that the rotor is not rotating.</i></p> <p>There is a plug on the right-hand side, which must be unscrewed counterclockwise using the emergency lid release key (catalog no. 18640). Then pull on the cap until the cover is open.</p> <p>The emergency opening of the cover can be used, for example, in the event of a power failure, failure of the control panel, etc.</p>
---	--

13. Guarantee


Manufacturer grants to the Buyer the guarantee on conditions specified in the Guarantee Certificate. Buyer forfeits the right to guarantee repair when using the device inconsistently with the User manual provisions, when damage results from the User's fault.

Repairs should be carried out in authorized service workshops, granted with the MPW Certificate.

The centrifuge shall be sent to repair after decontaminating disinfections. Information about authorized service workshops could be obtained from the Manufacturer.

	<ul style="list-style-type: none">▪ Guarantee period amounts to 24 months (unless otherwise specified in the purchase documents).▪ Guarantee conditions are described in guaranteed card.▪ The service life of the centrifuge specified by the manufacturer amounts to 10 years.▪ After 24 months from the start of the warranty period (date of purchase), a technical inspection of the centrifuge should be carried out (validation) by an authorized service of the manufacturer. Subsequent inspections should be carried out at annual intervals.▪ Maximum period of storage of not used centrifuge amounts to 1 year. After this period, a service authorized by manufacturer should carry out technical inspection of the centrifuge.▪ Manufacturer reserves the right to make technical changes in manufactured products.
---	---


14. Transport and storage

	<p>CAUTION! Due to the heavy weight of the device, lifting and carrying it may cause injury to the spine.</p>
<ul style="list-style-type: none">▪ Store the device only in a closed and dry room.▪ Remove rotor from centrifuge before transport.▪ Lift and carry with the adequate number of people.▪ Use transport equipment.▪ Use the original packaging and transport protection for transport.	

Transport and storage conditions.

	Storage (in the package)	Storage (without the package)	Transport
Temperature	-25 ÷ +55 °C	-5 ÷ +45 °C	-25 ÷ +60 °C (general) -20 ÷ +55 °C (air)
Relative humidity	10 ÷ 75 %	10 ÷ 75 %	10 ÷ 75 %
Pressure	70 ÷ 106 kPa	70 ÷ 106 kPa	30 ÷ 106 kPa

15. Disposal



- Dispose of the device in accordance with the applicable legal regulations in the country of use.
- In the countries of the European Community, the disposal of electrical equipment is regulated under the EU Directive 2012/19/EU on waste electrical and electronic equipment (WEEE). According to these regulations, centrifuges may not be collected together with municipal or household waste.
- Disposal regulations in individual EU countries may differ. In case of doubt, please contact the supplier of the device.

16. List of changes in the manual

<i>Rev.</i>	<i>Release date</i>	<i>Description of changes</i>
4	03.04.2023	Addition of markings used in the manual and on the device. Update of nameplate, CE declaration and equipment lists. Updating records regarding the intended use and disposal of the product. Removal of the RTC function.
5	16.06.2023	Removal of the USB communication function. Updating of the description in the technical data table. Updating the CE declaration of conformity, equipment list and nameplate.
6	11.01.2024	Updated "Installation" chapter, equipment list and name plate.

17. Manufacturer's info

"MPW MED. INSTRUMENTS" SPÓŁDZIELNIA PRACY

Boremłowska 46 Street

04-347 Warsaw

tel. (+48) 22 610 56 67 (sales department - POLAND)
(+48) 22 879 70 46 (sales department - outside POLAND)
(+48) 22 610 81 07 (service)
fax: (+48) 22 610 55 36
e-mail: mpw@mpw.pl
website: www.mpw.pl

000042924 - number of entries in the Waste Database

PL/CA01-01782 - identification number given by Office for Registration of Medicinal Products, Medical Devices and Biocidal Products.

Distributor's info

DISTRIBUTOR:

18. ANNEXES

A. Wyposażenie dodatkowe/Optional accessories

MPW-260/R/RH

WIRNIK / ROTOR

PARAMETRY/PARAMETERS (RCF [x g], Rmax [mm], α [°])

POJEMNIK/BUCKET

WKŁADKA / ADAPTER

[liczba probówek na wirnik/tubes per rotor] PROBÓWKA / TUBE

11199

RPM 18000 RCF 24270 Rmax 67 α 45

bez pojemnika/without bucket

bez wkładki/without adapter

[12] * 2-1,5 ml probówka (10,8x41,8 mm), Eppendorf®; [15011], 2 ml (10,8x41,8 mm); [15128], 1,5ml (10,8x40,5 mm)

14084

[12] * 0,5 ml probówka PCR (7,8 x 31 mm)
0,5 ml PCR tube (7,8 x 31 mm)

14126

[12] * 0,4 ml probówka PCR (5,7 x 48,6 mm)
0,4 ml PCR tube (5,7 x 48,6 mm)

14133

[12] * 0,2 ml probówka PCR (6 x 21,6 mm)
0,2 ml PCR tube (6 x 21,6 mm)

11213

RPM 5500 RCF 4227 Rmax 125 α 30

13276

bez wkładki/without adapter

[8] 15051 50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm)
50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm)

[8] * 50 ml probówka z dnem stożkowym z zakrętką (30 x 117 mm), Falcon®; [15052] 50ml (30 x 117mm)
50 ml tube, conical bottom, with cap (30 x 117 mm), Falcon®; [15052] 50ml Sarstedt® (30 x 117 mm)

[8] * 50 ml probówka z dnem stożkowym bez rantu (30 x 115 mm), Greiner®
50 ml tube, conical bottom, without skirt (30 x 115 mm), Greiner®

[8] * 50 ml probówka Advanced Oak Ridge (29x102 mm), Herolab® nr 25 32 11
50 ml tube, Advanced Oak Ridge (29 x 102 mm), Herolab® no. 25 32 11

14035

[8] 15046 14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt®
14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®

[8] 15053 10 ml probówka z pokrywką (16 x 106 mm)
10 ml tube with cap (16 x 106 mm)

[8] 15118 10 ml probówka szklana (16 x 100 mm)
10 ml glass tube (16 x 100 mm)

RCF max.=3000 RPM max.=4633

[8] * 15 ml Thermo Nalgene® (16 x 113 mm)
15 ml Thermo Nalgene® (16 x 113 mm)

14036

[8] * BD Vacutainer® (13 x 100 mm), (4-7 ml)

[8] * Greiner Vacuette® (13 x 100 mm), (3,5-6 ml)

[8] * 7 ml probówka szklana (12 x 100 mm)
7 ml glass tube (12 x 100 mm)

RCF max.=3000 RPM max.=4633

[8] * 6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt®
6 ml tube with cap (11,5 x 92 mm), Sarstedt®

14043

[8] * Greiner Vacuette® (13 x 75 mm), (1-4,5 ml)

[8] * Sarstedt S-Monovette® (13 x 75 mm), (2,7; 3; 4,3 ml)

[8] * Sarstedt S-Monovette® (13 x 90 mm), (4,9; 5,6 ml)

[8] * 5 ml probówka szklana (12 x 75 mm)
5 ml glass tube (12 x 75 mm)

RCF max.=3000 RPM max.=4633

[8] * 5 ml probówka z korkiem (12 x 85 mm), Sarstedt®
5 ml tube with cap (12 x 85 mm), Sarstedt®

14071

[8] 15055 30 ml probówka z pokrywką (25,4 x 103,2 mm)
30 ml tube with cap (25,4 x 103,2 mm)

[8] * 28 ml Thermo Nalgene® Oak Ridge (25,4 x 101,8 mm)

[8] * 30 ml probówka z pokrywką (25,5 x 94 mm), Nalgene®
30 ml tube with cap (25,5 x 94 mm), Nalgene®

[8] * 30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm)
30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm)

* probówka niedostępna w ofercie MPW lub dostępny odpowiednik (np.[15050]), patrz kolumna z prawej
tube is not offered by MPW or equivalent is available (e.g. [15050]), see column on the right

A. Wyposażenie dodatkowe/Optional accessories
MPW-260/R/RH
14073

[8] 15046 14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt®
14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®

[8] 15053 10 ml probówka z pokrywką (16 x 106 mm)
10 ml tube with cap (16 x 106 mm)

[8] * BD Vacutainer® (16 x 100 mm), (2,5-11 ml)

[8] * Greiner Vacuette® (16 x 100 mm), (7-9 ml)

[8] * Sarstedt S-Monovette® (15 x 92 mm), (7,5; 8,2; 8,5 ml)

[8] * Sarstedt S-Monovette® (16 x 92 mm), (9; 10 ml)

[8] 15118 10 ml probówka szklana (16 x 100 mm)
10 ml glass tube (16 x 100 mm)
RCF max.=3000 RPM max.=4633

14089

[8] * 15 ml probówka z dnem stożkowym z zakrętką (17 x 120 mm), Falcon®; [15050], 15ml (17 x 120 mm)
15 ml tube, conical bottom, with cap (17 x 120 mm), Falcon®; [15050] 15ml Sarstedt®(17 x 120 mm)

14248

[8] 15055 30 ml probówka z pokrywką (25,4 x 103,2 mm)
30 ml tube with cap (25,4 x 103,2 mm)

14089+14868

[8] * 5 ml probówka z korkiem wciskany (17 x 54,2 mm), Eppendorf®
5 ml tube with snap cap (17 x 54,2 mm), Eppendorf®

[8] * 5 ml probówka z korkiem zakręcany (17 x 66 mm), Eppendorf®
5 ml tube with screw cap (17 x 66 mm), Eppendorf®

13278+17151
bez wkładki/without adapter

[8] 15051 50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm)
50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm)

[8] * 50 ml probówka z dnem stożkowym z zakrętką (30 x 117 mm), Falcon®; [15052] 50ml (30 x 117mm)
50 ml tube, conical bottom, with cap (30 x 117 mm), Falcon®; [15052] 50ml Sarstedt® (30 x 117 mm)

[8] * 50 ml probówka z dnem stożkowym bez rantu (30 x 115 mm), Greiner®
50 ml tube, conical bottom, without skirt (30 x 115 mm), Greiner®

[8] * 50 ml probówka Advanced Oak Ridge (29x102 mm), Herolab® nr 25 32 11
50 ml tube, Advanced Oak Ridge (29 x 102 mm), Herolab® no. 25 32 11

14035

[8] 15046 14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt®
14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®

[8] 15053 10 ml probówka z pokrywką (16 x 106 mm)
10 ml tube with cap (16 x 106 mm)

[8] 15118 10 ml probówka szklana (16 x 100 mm)
10 ml glass tube (16 x 100 mm)
RCF max.=3000 RPM max.=4633

[8] * 15 ml Thermo Nalgene® (16 x 113 mm)
15 ml Thermo Nalgene® (16 x 113 mm)

14036

[8] * 7 ml probówka szklana (12 x 100 mm)
7 ml glass tube (12 x 100 mm)
RCF max.=3000 RPM max.=4633

[8] * 6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt®
6 ml tube with cap (11,5 x 92 mm), Sarstedt®

14043

[8] * Greiner Vacuette® (13 x 75 mm), (1-4,5 ml)

[8] * 5 ml probówka szklana (12 x 75 mm)
5 ml glass tube (12 x 75 mm)
RCF max.=3000 RPM max.=4633

[8] * 5 ml probówka z korkiem (12 x 85 mm), Sarstedt®
5 ml tube with cap (12 x 85 mm), Sarstedt®

14071

[8] 15055 30 ml probówka z pokrywką (25,4 x 103,2 mm)
30 ml tube with cap (25,4 x 103,2 mm)

[8] * 28 ml Thermo Nalgene® Oak Ridge (25,4 x 101,8 mm)

[8] * 30 ml probówka z pokrywką (25,5 x 94 mm), Nalgene®
30 ml tube with cap (25,5 x 94 mm), Nalgene®

[8] * 30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm)
30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm)

14073

[8] 15046 14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt®
14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®

[8] 15053 10 ml probówka z pokrywką (16 x 106 mm)
10 ml tube with cap (16 x 106 mm)

[8] * BD Vacutainer® (16 x 100 mm), (2,5-11 ml)

[8] * Greiner Vacuette® (16 x 100 mm), (7-9 ml)

[8] * Sarstedt S-Monovette® (15 x 92 mm), (7,5; 8,2; 8,5 ml)

[8] * Sarstedt S-Monovette® (16 x 92 mm), (9; 10 ml)

[8] 15118 10 ml probówka szklana (16 x 100 mm)
10 ml glass tube (16 x 100 mm)
RCF max.=3000 RPM max.=4633

* probówka niedostępna w ofercie MPW lub dostępny odpowiednik (np.[15050]), patrz kolumna z prawej
tube is not offered by MPW or equivalent is available (e.g. [15050]), see column on the right

A. Wyposażenie dodatkowe/Optional accessories
MPW-260/R/RH
14089

- [8] * 15 ml probówka z dnem stożkowym z zakrętką (17 x 120 mm), Falcon®; [15050], 15ml (17 x 120 mm)
15 ml tube, conical bottom, with cap (17 x 120 mm), Falcon®; [15050] 15ml Sarstedt®(17 x 120 mm)

14248

- [8] 15055 30 ml probówka z pokrywką (25,4 x 103,2 mm)
30 ml tube with cap (25,4 x 103,2 mm)

14089+14868

- [8] * 5 ml probówka z korkiem wciskany (17 x 54,2 mm), Eppendorf®
5 ml tube with snap cap (17 x 54,2 mm), Eppendorf®

11216

RPM 14000 RCF 19064 Rmax 87 \pm 45

bez pojemnika/without bucket

bez wkładki/without adapter

- [12] * 5 ml probówka z korkiem (12 x 85 mm), Sarstedt®
5 ml tube with cap (12 x 85 mm), Sarstedt®

11217

RPM 6000 RCF 4226 Rmax 105 \pm 30

13080
14082

- [10] * BD Vacutainer® (13 x 100 mm), (4-7 ml)
[10] * Greiner Vacuette® (13 x 100 mm), (3,5-6 ml)
[10] * Sarstedt S-Monovette® (11 x 92 mm), (4,5; 5 ml)
[10] * 7 ml probówka szklana (12 x 100 mm)
7 ml glass tube (12 x 100 mm)
RCF max.=3000 RPM max.=5055

- [10] * 6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt®
6 ml tube with cap (11,5 x 92 mm), Sarstedt®

bez wkładki/without adapter

- [10] 15046 14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt®
14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®

- [10] 15053 10 ml probówka z pokrywką (16 x 106 mm)
10 ml tube with cap (16 x 106 mm)

- [10] * Greiner Vacuette® (16 x 100 mm), (7-9 ml)
[10] * Sarstedt S-Monovette® (15 x 92 mm), (7,5; 8,2; 8,5 ml)

- [10] * 15 ml probówka z dnem stożkowym z zakrętką (17 x 120 mm), Falcon®; [15050], 15ml (17 x 120 mm)
15 ml tube, conical bottom, with cap (17 x 120 mm), Falcon®; [15050] 15ml Sarstedt®(17 x 120 mm)

- [10] * BD Vacutainer® (16 x 100 mm), (2,5-11 ml)

- [10] * Sarstedt S-Monovette® (16 x 92 mm), (9; 10 ml)

- [10] 15118 10 ml probówka szklana (16 x 100 mm)
10 ml glass tube (16 x 100 mm)

RCF max.=3000 RPM max.=5055

- [10] * 15 ml Thermo Nalgene® (16 x 113 mm)
15 ml Thermo Nalgene® (16 x 113 mm)

14082+14815

- [10] * 5 ml probówka szklana (12 x 75 mm)
5 ml glass tube (12 x 75 mm)
RCF max.=3000 RPM max.=5554

14082+14815 Rmax 87 RCF 3502

- [10] * Greiner Vacuette® (13 x 75 mm), (1-4,5 ml)
[10] * Sarstedt S-Monovette® (13 x 75 mm), (2,7; 3; 4,3 ml)
[10] * BD Vacutainer® (13 x 75 mm), (1,6-5,3 ml)
[10] * Sarstedt S-Monovette® (11 x 66 mm), (1,6; 2; 2,7; 3; 3,1 ml)
[10] * Sarstedt S-Monovette® (13 x 65 mm), (2,6; 2,9; 3,4; 3,8 ml)

14815 Rmax 87 RCF 3502

- [10] 15121 10 ml probówka z dnem okrągłym i pokrywką (17 x 70 mm)
10 ml tube, round bottom, with cap (17 x 70 mm)

- [10] * Sarstedt S-Monovette® (15 x 75 mm), (4; 4,3; 5,5 ml)

- [10] * 10 ml Thermo Nalgene® Oak Ridge (16 x 81,5 mm)

11461

RPM 15100 RCF 21158 Rmax 83 \pm 45

bez pojemnika/without bucket

bez wkładki/without adapter

- [24] * 2-1,5 ml probówka (10,8x41,8 mm), Eppendorf®; [15011], 2 ml (10,8x41,8 mm); [15128], 1,5ml (10,8x40,5 mm)

14084

- [24] * 0,5 ml probówka PCR (7,8 x 31 mm)
0,5 ml PCR tube (7,8 x 31 mm)

14126

- [24] * 0,4 ml probówka PCR (5,7 x 48,6 mm)
0,4 ml PCR tube (5,7 x 48,6 mm)

14133

- [24] * 0,2 ml probówka PCR (6 x 21,6 mm)
0,2 ml PCR tube (6 x 21,6 mm)

A. Wyposażenie dodatkowe/Optional accessories
MPW-260/R/RH
11462
RPM 14000 RCF 18188 Rmax 83 \pm 45
bez pojemnika/without bucket
bez wkładki/without adapter

[36] * 2-1,5 ml probówka (10,8x41,8 mm), Eppendorf®; [15011], 2 ml (10,8x41,8 mm); [15128], 1,5ml (10,8x40,5 mm)

14084

 [36] * 0,5 ml probówka PCR (7,8 x 31 mm)
 0,5 ml PCR tube (7,8 x 31 mm)

14126

 [36] * 0,4 ml probówka PCR (5,7 x 48,6 mm)
 0,4 ml PCR tube (5,7 x 48,6 mm)

14133

 [36] * 0,2 ml probówka PCR (6 x 21,6 mm)
 0,2 ml PCR tube (6 x 21,6 mm)

11501
RPM 4500 RCF 2966 Rmax 131 \pm 30
13080
14082

[30] * BD Vacutainer® (13 x 100 mm), (4-7 ml)

[30] * Greiner Vacuette® (13 x 100 mm), (3,5-6 ml)

[30] * Sarstedt S-Monovette® (11 x 92 mm), (4,5; 5 ml)

 [30] * 7 ml probówka szklana (12 x 100 mm)
 7 ml glass tube (12 x 100 mm)

 [30] * 6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt®
 6 ml tube with cap (11,5 x 92 mm), Sarstedt®

bez wkładki/without adapter

 [30] 15046 14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt®
 14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®

 [30] 15053 10 ml probówka z pokrywką (16 x 106 mm)
 10 ml tube with cap (16 x 106 mm)

 [30] * 15 ml probówka z dnem stożkowym z zakrętką (17 x 120 mm), Falcon®; [15050], 15ml (17 x 120 mm)
 15 ml tube, conical bottom, with cap (17 x 120 mm), Falcon®; [15050] 15ml Sarstedt®(17 x 120 mm)

[30] * BD Vacutainer® (16 x 100 mm), (2,5-11 ml)

[30] * Greiner Vacuette® (16 x 100 mm), (7-9 ml)

[30] * Sarstedt S-Monovette® (15 x 92 mm), (7,5; 8,2; 8,5 ml)

[30] * Sarstedt S-Monovette® (16 x 92 mm), (9; 10 ml)

 [30] 15118 10 ml probówka szklana (16 x 100 mm)
 10 ml glass tube (16 x 100 mm)

 [30] * 15 ml Thermo Nalgene® (16 x 113 mm)
 15 ml Thermo Nalgene® (16 x 113 mm)

14082+14815 Rmax 120 RCF 2717

[30] * BD Vacutainer® (13 x 75 mm), (1,6-5,3 ml)

[30] * Greiner Vacuette® (13 x 75 mm), (1-4,5 ml)

[30] * Sarstedt S-Monovette® (11 x 66 mm), (1,6; 2; 2,7; 3; 3,1 ml)

[30] * Sarstedt S-Monovette® (13 x 65 mm), (2,6; 2,9; 3,4; 3,8 ml)

[30] * Sarstedt S-Monovette® (13 x 75 mm), (2,7; 3; 4,3 ml)

 [30] * 5 ml probówka szklana (12 x 75 mm)
 5 ml glass tube (12 x 75 mm)

14815 Rmax 120 RCF 2717

 [30] 15121 10 ml probówka z dnem okrągłym i pokrywką (17 x 70 mm)
 10 ml tube, round bottom, with cap (17 x 70 mm)

[30] * Sarstedt S-Monovette® (15 x 75 mm), (4; 4,3; 5,5 ml)

[30] * 10 ml Thermo Nalgene® Oak Ridge (16 x 81,5 mm)

11715
RPM 14000 RCF 15558 Rmax 71 \pm 30
bez pojemnika/without bucket
bez wkładki/without adapter

 [10] 15121 10 ml probówka z dnem okrągłym i pokrywką (17 x 70 mm)
 10 ml tube, round bottom, with cap (17 x 70 mm)

11716
RPM 14000 RCF 15339 Rmax 70 \pm 45
bez pojemnika/without bucket
bez wkładki/without adapter

 [4] * 8 x 0,2 ml probówki szeregowe PCR-strip (10,2 x 72,4 mm)
 8 x 0,2 ml PCR strip (10,2 x 72,4 mm)

 [32] * 0,2 ml probówka PCR (6 x 21,6 mm)
 0,2 ml PCR tube (6 x 21,6 mm)

 [4] * 8 x 0,2 ml probówki szeregowe PCR strip (7,3 x 77,2 mm)
 8 x 0,2 ml PCR strip (7,3 x 77,2 mm)

 [4] * 4 x 0,2 ml probówki szeregowe PCR-strip (10,2 x 37,2 mm)
 4 x 0,2 ml PCR strip (10,2 x 37,2 mm)

 * probówka niedostępna w ofercie MPW lub dostępny odpowiednik (np.[15050]), patrz kolumna z prawej
 tube is not offered by MPW or equivalent is available (e.g. [15050]), see column on the right

A. Wyposażenie dodatkowe/Optional accessories
MPW-260/R/RH
11718
RPM 6300 RCF 5014 Rmax 113 \pm 30
13719
14024

[4] * 15 ml probówka z dnem stożkowym z zakrętką (17 x 120 mm), Falcon®; [15050], 15ml (17 x 120 mm)
15 ml tube, conical bottom, with cap (17 x 120 mm), Falcon®; [15050] 15ml Sarstedt®(17 x 120 mm)

14196

[4] 15040 100 ml probówka z pokrywką (45,2 x 103,7 mm)
100 ml tube with cap (45,2 x 103,7 mm)

14224

[4] 15055 30 ml probówka z pokrywką (25,4 x 103,2 mm)
30 ml tube with cap (25,4 x 103,2 mm)

[4] 15222 30 ml probówka z pokrywką (25 x 94mm), Sterilin®
30 ml tube with cap (25 x 94 mm), Sterilin®

[4] 15223 30 ml probówka z pokrywką (25 x 94 mm), Sterilin®
30 ml tube with cap (25 x 94 mm), Sterilin®

[4] * 30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm)
30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm)

14226

[4] * 50 ml probówka z dnem stożkowym z rantem (30 x 115 mm), Greiner®
50 ml tube, conical bottom, skirted (30 x 115 mm), Greiner®

14189+14188

[4] 15051 50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm)
50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm)

[4] * 50 ml probówka z dnem stożkowym z zakrętką (30 x 117 mm), Falcon®; [15052] 50ml (30 x 117mm)
50 ml tube, conical bottom, with cap (30 x 117 mm), Falcon®; [15052] 50ml Sarstedt® (30 x 117 mm)

[4] * 50 ml probówka z dnem stożkowym bez rantu (30 x 115 mm), Greiner®
50 ml tube, conical bottom, without skirt (30 x 115 mm), Greiner®

[4] * 50 ml probówka Advanced Oak Ridge (29x102 mm), Herolab® nr 25 32 11
50 ml tube, Advanced Oak Ridge (29 x 102 mm), Herolab® no. 25 32 11

14190+14188

[4] 15055 30 ml probówka z pokrywką (25,4 x 103,2 mm)
30 ml tube with cap (25,4 x 103,2 mm)

11740
RPM 5500 RCF 4058 Rmax 120 \pm 30
13080
14082

[12] * BD Vacutainer® (13 x 100 mm), (4-7 ml)

[12] * Greiner Vacuette® (13 x 100 mm), (3,5-6 ml)

[12] * Sarstedt S-Monovette® (11 x 92 mm), (4,5; 5 ml)

[12] * 7 ml probówka szklana (12 x 100 mm)
7 ml glass tube (12 x 100 mm)
RCF max.=3000 RPM max.=4729

[12] * 6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt®
6 ml tube with cap (11,5 x 92 mm), Sarstedt®

bez wkładki/without adapter

[12] 15046 14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt®
14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®

[12] 15053 10 ml probówka z pokrywką (16 x 106 mm)
10 ml tube with cap (16 x 106 mm)

[12] * 15 ml probówka z dnem stożkowym z zakrętką (17 x 120 mm), Falcon®; [15050], 15ml (17 x 120 mm)
15 ml tube, conical bottom, with cap (17 x 120 mm), Falcon®; [15050] 15ml Sarstedt®(17 x 120 mm)

[12] * BD Vacutainer® (16 x 100 mm), (2,5-11 ml)

[12] * Greiner Vacuette® (16 x 100 mm), (7-9 ml)

[12] * Sarstedt S-Monovette® (15 x 92 mm), (7,5; 8,2; 8,5 ml)

[12] * Sarstedt S-Monovette® (16 x 92 mm), (9; 10 ml)

[12] 15118 10 ml probówka szklana (16 x 100 mm)
10 ml glass tube (16 x 100 mm)
RCF max.=3000 RPM max.=4729

[12] * 15 ml Thermo Nalgene® (16 x 113 mm)

15 ml Thermo Nalgene® (16 x 113 mm)

14082+14815

[12] * 5 ml probówka szklana (12 x 75 mm)
5 ml glass tube (12 x 75 mm)
RCF max.=3000 RPM max.=5154

14082+14815 Rmax 101 RCF 3416

[12] * BD Vacutainer® (13 x 75 mm), (1,6-5,3 ml)

[12] * Greiner Vacuette® (13 x 75 mm), (1-4,5 ml)

[12] * Sarstedt S-Monovette® (11 x 66 mm), (1,6; 2; 2,7; 3; 3,1 ml)

[12] * Sarstedt S-Monovette® (13 x 65 mm), (2,6; 2,9; 3,4; 3,8 ml)

[12] * Sarstedt S-Monovette® (13 x 75 mm), (2,7; 3; 4,3 ml)

14815 Rmax 101 RCF 3416

[12] 15121 10 ml probówka z dnem okrągłym i pokrywką (17 x 70 mm)

10 ml tube, round bottom, with cap (17 x 70 mm)

[12] * Sarstedt S-Monovette® (15 x 75 mm), (4; 4,3; 5,5 ml)

* probówka niedostępna w ofercie MPW lub dostępny odpowiednik (np.[15050]), patrz kolumna z prawej
tube is not offered by MPW or equivalent is available (e.g. [15050]), see column on the right

A. Wyposażenie dodatkowe/Optional accessories
MPW-260/R/RH

[12] * 10 ml Thermo Nalgene® Oak Ridge (16 x 81,5 mm)

11743
RPM 4500 RCF 2604 Rmax 115 α 30
13329
bez wkładki/without adapter

 [12] 15055 30 ml probówka z pokrywką (25,4 x 103,2 mm)
 30 ml tube with cap (25,4 x 103,2 mm)

 [12] 15222 30 ml probówka z pokrywką (25 x 94mm), Sterilin®
 30 ml tube with cap (25 x 94 mm), Sterilin®

 [12] 15223 30 ml probówka z pokrywką (25 x 94 mm), Sterilin®
 30 ml tube with cap (25 x 94 mm), Sterilin®

 [12] * 30 ml probówka z pokrywką (25,5 x 94 mm), Nalgene®
 30 ml tube with cap (25,5 x 94 mm), Nalgene®

14256

 [12] 15046 14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt®
 14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®

 [12] 15053 10 ml probówka z pokrywką (16 x 106 mm)
 10 ml tube with cap (16 x 106 mm)

 [12] 15118 10 ml probówka szklana (16 x 100 mm)
 10 ml glass tube (16 x 100 mm)

 [12] * 15 ml Thermo Nalgene® (16 x 113 mm)
 15 ml Thermo Nalgene® (16 x 113 mm)

14255

[12] * Sarstedt S-Monovette® (11 x 92 mm), (4,5; 5 ml)

 [12] * 7 ml probówka szklana (12 x 100 mm)
 7 ml glass tube (12 x 100 mm)

11744
RPM 4500 RCF 2830 Rmax 125 α 30
13276
bez wkładki/without adapter

 [10] 15051 50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm)
 50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm)

 [10] * 50 ml probówka z dnem stożkowym z zakrętką (30 x 117 mm), Falcon®; [15052] 50ml (30 x 117mm)
 50 ml tube, conical bottom, with cap (30 x 117 mm), Falcon®; [15052] 50ml Sarstedt® (30 x 117 mm)

 [10] * 50 ml probówka z dnem stożkowym bez rantu (30 x 115 mm), Greiner®
 50 ml tube, conical bottom, without skirt (30 x 115 mm), Greiner®

 [10] * 50 ml probówka Advanced Oak Ridge (29x102 mm), Herolab® nr 25 32 11
 50 ml tube, Advanced Oak Ridge (29 x 102 mm), Herolab® no. 25 32 11

14035

 [10] 15046 14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt®
 14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®

 [10] 15053 10 ml probówka z pokrywką (16 x 106 mm)
 10 ml tube with cap (16 x 106 mm)

 [10] 15118 10 ml probówka szklana (16 x 100 mm)
 10 ml glass tube (16 x 100 mm)

 [10] * 15 ml Thermo Nalgene® (16 x 113 mm)
 15 ml Thermo Nalgene® (16 x 113 mm)

14036

[10] * BD Vacutainer® (13 x 100 mm), (4-7 ml)

[10] * Greiner Vacuette® (13 x 100 mm), (3,5-6 ml)

 [10] * 7 ml probówka szklana (12 x 100 mm)
 7 ml glass tube (12 x 100 mm)

 [10] * 6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt®
 6 ml tube with cap (11,5 x 92 mm), Sarstedt®

14043

[10] * Greiner Vacuette® (13 x 75 mm), (1-4,5 ml)

[10] * Sarstedt S-Monovette® (13 x 75 mm), (2,7; 3; 4,3 ml)

[10] * Sarstedt S-Monovette® (13 x 90 mm), (4,9; 5,6 ml)

 [10] * 5 ml probówka z korkiem (12 x 85 mm), Sarstedt®
 5 ml tube with cap (12 x 85 mm), Sarstedt®

 [10] * 5 ml probówka szklana (12 x 75 mm)
 5 ml glass tube (12 x 75 mm)

14071

 [10] 15055 30 ml probówka z pokrywką (25,4 x 103,2 mm)
 30 ml tube with cap (25,4 x 103,2 mm)

[10] * 28 ml Thermo Nalgene® Oak Ridge (25,4 x 101,8 mm)

 [10] * 30 ml probówka z pokrywką (25,5 x 94 mm), Nalgene®
 30 ml tube with cap (25,5 x 94 mm), Nalgene®

 [10] * 30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm)
 30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm)

14073

 [10] 15046 14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt®
 14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®

 [10] 15053 10 ml probówka z pokrywką (16 x 106 mm)
 10 ml tube with cap (16 x 106 mm)

* probówka niedostępna w ofercie MPW lub dostępny odpowiednik (np.[15050]), patrz kolumna z prawej

tube is not offered by MPW or equivalent is available (e.g. [15050]), see column on the right

A. Wyposażenie dodatkowe/Optional accessories
MPW-260/R/RH

[10]	*	Greiner Vacuette® (16 x 100 mm), (7-9 ml)
[10]	*	Sarstedt S-Monovette® (15 x 92 mm), (7,5; 8,2; 8,5 ml)
[10]	*	BD Vacutainer® (16 x 100 mm), (2,5-11 ml)
[10]	*	Sarstedt S-Monovette® (16 x 92 mm), (9; 10 ml)
[10]	15118	10 ml probówka szklana (16 x 100 mm) 10 ml glass tube (16 x 100 mm)
[10]	*	15 ml Thermo Nalgene® (16 x 113 mm) 15 ml Thermo Nalgene® (16 x 113 mm)
14089		
[10]	*	15 ml probówka z dnem stożkowym z zakrętką (17 x 120 mm), Falcon®; [15050], 15ml (17 x 120 mm) 15 ml tube, conical bottom, with cap (17 x 120 mm), Falcon®; [15050] 15ml Sarstedt®(17 x 120 mm)
14248		
[10]	15055	30 ml probówka z pokrywką (25,4 x 103,2 mm) 30 ml tube with cap (25,4 x 103,2 mm)
14089+14868		
[10]	*	5 ml probówka z korkiem wciskany (17 x 54,2 mm), Eppendorf® 5 ml tube with snap cap (17 x 54,2 mm), Eppendorf®
[10]	*	5 ml probówka z korkiem zakręcany (17 x 66 mm), Eppendorf® 5 ml tube with screw cap (17 x 66 mm), Eppendorf®
11745		
RPM 5000 RCF 3354 Rmax 120 \pm 30		
13080		
14082		
[24]	*	BD Vacutainer® (13 x 100 mm), (4-7 ml)
[24]	*	Greiner Vacuette® (13 x 100 mm), (3,5-6 ml)
[24]	*	Sarstedt S-Monovette® (11 x 92 mm), (4,5; 5 ml)
[24]	*	7 ml probówka szklana (12 x 100 mm) 7 ml glass tube (12 x 100 mm) RCF max.=3000 RPM max.=4729
[24]	*	6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt® 6 ml tube with cap (11,5 x 92 mm), Sarstedt® bez wkładki/without adapter
[24]	15046	14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt® 14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®
[24]	15053	10 ml probówka z pokrywką (16 x 106 mm) 10 ml tube with cap (16 x 106 mm)
[24]	*	15 ml probówka z dnem stożkowym z zakrętką (17 x 120 mm), Falcon®; [15050], 15ml (17 x 120 mm) 15 ml tube, conical bottom, with cap (17 x 120 mm), Falcon®; [15050] 15ml Sarstedt®(17 x 120 mm)
[24]	*	BD Vacutainer® (16 x 100 mm), (2,5-11 ml)
[24]	*	Greiner Vacuette® (16 x 100 mm), (7-9 ml)
[24]	*	Sarstedt S-Monovette® (15 x 92 mm), (7,5; 8,2; 8,5 ml)
[24]	*	Sarstedt S-Monovette® (16 x 92 mm), (9; 10 ml)
[24]	15118	10 ml probówka szklana (16 x 100 mm) 10 ml glass tube (16 x 100 mm) RCF max.=3000 RPM max.=4729
[24]	*	15 ml Thermo Nalgene® (16 x 113 mm) 15 ml Thermo Nalgene® (16 x 113 mm)
14082+14815 Rmax 105 RCF 2935		
[24]	*	BD Vacutainer® (13 x 75 mm), (1,6-5,3 ml)
[24]	*	Greiner Vacuette® (13 x 75 mm), (1-4,5 ml)
[24]	*	Sarstedt S-Monovette® (11 x 66 mm), (1,6; 2; 2,7; 3; 3,1 ml)
[24]	*	Sarstedt S-Monovette® (13 x 65 mm), (2,6; 2,9; 3,4; 3,8 ml)
[24]	*	Sarstedt S-Monovette® (13 x 75 mm), (2,7; 3; 4,3 ml)
[24]	*	5 ml probówka szklana (12 x 75 mm) 5 ml glass tube (12 x 75 mm)
14815 Rmax 105 RCF 2935		
[24]	15121	10 ml probówka z dnem okrągłym i pokrywką (17 x 70 mm) 10 ml tube, round bottom, with cap (17 x 70 mm)
[24]	*	Sarstedt S-Monovette® (15 x 75 mm), (4; 4,3; 5,5 ml)
[24]	*	10 ml Thermo Nalgene® Oak Ridge (16 x 81,5 mm)
11746		
RPM 6000 RCF 4427 Rmax 110 \pm 30		
13276		
bez wkładki/without adapter		
[6]	*	50 ml probówka z dnem stożkowym z zakrętką (30 x 117 mm), Falcon®; [15052] 50ml (30 x 117mm) 50 ml tube, conical bottom, with cap (30 x 117 mm), Falcon®; [15052] 50ml Sarstedt® (30 x 117 mm)
[6]	*	50 ml probówka z dnem stożkowym bez rantu (30 x 115 mm), Greiner® 50 ml tube, conical bottom, without skirt (30 x 115 mm), Greiner®
[6]	15051	50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm) 50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm)
[6]	*	50 ml probówka Advanced Oak Ridge (29x102 mm), Herolab® nr 25 32 11 50 ml tube, Advanced Oak Ridge (29 x 102 mm), Herolab® no. 25 32 11
14035		
[6]	15046	14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt® 14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®

* probówka niedostępna w ofercie MPW lub dostępny odpowiednik (np.[15050]), patrz kolumna z prawej
tube is not offered by MPW or equivalent is available (e.g. [15050]), see column on the right

A. Wyposażenie dodatkowe/Optional accessories
MPW-260/R/RH

[6]	15053	10 ml probówka z pokrywką (16 x 106 mm) 10 ml tube with cap (16 x 106 mm)
[6]	15118	10 ml probówka szklana (16 x 100 mm) 10 ml glass tube (16 x 100 mm) RCF max.=3000 RPM max.=4939
[6]	*	15 ml Thermo Nalgene® (16 x 113 mm) 15 ml Thermo Nalgene® (16 x 113 mm)
14036		
[6]	*	BD Vacutainer® (13 x 100 mm), (4-7 ml)
[6]	*	Greiner Vacuette® (13 x 100 mm), (3,5-6 ml)
[6]	*	7 ml probówka szklana (12 x 100 mm) 7 ml glass tube (12 x 100 mm) RCF max.=3000 RPM max.=4939
[6]	*	6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt® 6 ml tube with cap (11,5 x 92 mm), Sarstedt®
14043		
[6]	*	Greiner Vacuette® (13 x 75 mm), (1-4,5 ml)
[6]	*	Sarstedt S-Monovette® (13 x 75 mm), (2,7; 3; 4,3 ml)
[6]	*	Sarstedt S-Monovette® (13 x 90 mm), (4,9; 5,6 ml)
[6]	*	5 ml probówka z korkiem (12 x 85 mm), Sarstedt® 5 ml tube with cap (12 x 85 mm), Sarstedt®
[6]	*	5 ml probówka szklana (12 x 75 mm) 5 ml glass tube (12 x 75 mm) RCF max.=3000 RPM max.=4939
14071		
[6]	15055	30 ml probówka z pokrywką (25,4 x 103,2 mm) 30 ml tube with cap (25,4 x 103,2 mm)
[6]	*	28 ml Thermo Nalgene® Oak Ridge (25,4 x 101,8 mm)
[6]	*	30 ml probówka z pokrywką (25,5 x 94 mm), Nalgene® 30 ml tube with cap (25,5 x 94 mm), Nalgene®
[6]	*	30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm) 30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm)
14073		
[6]	15046	14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt® 14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®
[6]	15053	10 ml probówka z pokrywką (16 x 106 mm) 10 ml tube with cap (16 x 106 mm)
[6]	*	BD Vacutainer® (16 x 100 mm), (2,5-11 ml)
[6]	*	Greiner Vacuette® (16 x 100 mm), (7-9 ml)
[6]	*	Sarstedt S-Monovette® (15 x 92 mm), (7,5; 8,2; 8,5 ml)
[6]	*	Sarstedt S-Monovette® (16 x 92 mm), (9; 10 ml)
[6]	15118	10 ml probówka szklana (16 x 100 mm) 10 ml glass tube (16 x 100 mm) RCF max.=3000 RPM max.=4939
14089		
[6]	*	15 ml probówka z dnem stożkowym z zakrętką (17 x 120 mm), Falcon®; [15050], 15ml (17 x 120 mm) 15 ml tube, conical bottom, with cap (17 x 120 mm), Falcon®; [15050] 15ml Sarstedt®(17 x 120 mm)
14248		
[6]	15055	30 ml probówka z pokrywką (25,4 x 103,2 mm) 30 ml tube with cap (25,4 x 103,2 mm)
14089+14868		
[6]	*	5 ml probówka z korkiem wciskany (17 x 54,2 mm), Eppendorf® 5 ml tube with snap cap (17 x 54,2 mm), Eppendorf®
[6]	*	5 ml probówka z korkiem zakręcany (17 x 66 mm), Eppendorf® 5 ml tube with screw cap (17 x 66 mm), Eppendorf®
11760		
RPM 14600 RCF 20257 Rmax 85 ɳ 45		
bez pojemnika/without bucket		
bez wkładki/without adapter		
[24]	*	2-1,5 ml probówka (10,8x41,8 mm), Eppendorf®; [15011], 2 ml (10,8x41,8 mm); [15128], 1,5ml (10,8x40,5 mm)
[24]	*	2 ml probówki z filtrem - spin columns (10,8 x 46 mm) 2 ml spin columns (with filter) (10,8 x 46 mm); [15011], 2 ml (10,8x41,8 mm); [15128], 1,5ml
14084		
[24]	*	0,5 ml probówka PCR (7,8 x 31 mm) 0,5 ml PCR tube (7,8 x 31 mm)
14126		
[24]	*	0,4 ml probówka PCR (5,7 x 48,6 mm) 0,4 ml PCR tube (5,7 x 48,6 mm)
14133		
[24]	*	0,2 ml probówka PCR (6 x 21,6 mm) 0,2 ml PCR tube (6 x 21,6 mm)
11943		
RPM 12000 RCF 13684 Rmax 85 ɳ 45		
bez pojemnika/without bucket		

* probówka niedostępna w ofercie MPW lub dostępny odpowiednik (np.[15050]), patrz kolumna z prawej
tube is not offered by MPW or equivalent is available (e.g. [15050]), see column on the right

A. Wyposażenie dodatkowe/Optional accessories
MPW-260/R/RH
bez wkładki/without adapter

[20]	*	1,6 ml probówka Cryo (12,3 x 46,5 mm) 1,6 ml Cryo tube (12,3 x 46,5 mm)
[20]	*	1,8 ml probówka Cryo (12,3 x 46,5 mm) 1,8 ml Cryo tube (12,3 x 46,5 mm)

11944
RPM 12000 RCF 13684 Rmax 85 ϕ 45
bez pojemnika/without bucket
bez wkładki/without adapter

[6]	*	5 ml probówka z korkiem zakręcany (17 x 66 mm), Eppendorf® 5 ml tube with screw cap (17 x 66 mm), Eppendorf®
[12]	*	5 ml probówka z korkiem wciskany (17 x 54,2 mm), Eppendorf® 5 ml tube with snap cap (17 x 54,2 mm), Eppendorf®

12200
RPM 4000 RCF 2504 Rmax 140 ϕ 90
13200
14013

[32]	*	BD Vacutainer® (13 x 75 mm), (1,6-5,3 ml)
[32]	*	BD Vacutainer® (13 x 100 mm), (4-7 ml)
[32]	*	Greiner Vacuette® (13 x 75 mm), (1-4,5 ml)
[32]	*	Greiner Vacuette® (13 x 100 mm), (3,5-6 ml)
[32]	*	Sarstedt S-Monovette® (13 x 65 mm), (2,6; 2,9; 3,4; 3,8 ml)
[32]	*	Sarstedt S-Monovette® (13 x 75 mm), (2,7; 3; 4,3 ml)
[32]	*	Sarstedt S-Monovette® (11 x 92 mm), (4,5; 5 ml)
[32]	*	Sarstedt S-Monovette® (13 x 90 mm), (4,9; 5,6 ml)
[32]	*	7 ml probówka szklana (12 x 100 mm) 7 ml glass tube (12 x 100 mm)
[32]	*	5 ml probówka szklana (12 x 75 mm) 5 ml glass tube (12 x 75 mm)
[32]	*	5 ml probówka z korkiem (12 x 85 mm), Sarstedt® 5 ml tube with cap (12 x 85 mm), Sarstedt®
[32]	*	6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt® 6 ml tube with cap (11,5 x 92 mm), Sarstedt®

14016

[28]	15053	10 ml probówka z pokrywką (16 x 106 mm) 10 ml tube with cap (16 x 106 mm)
[28]	*	BD Vacutainer® (16 x 100 mm), (2,5-11 ml)
[28]	*	Greiner Vacuette® (16 x 100 mm), (7-9 ml)
[28]	*	Sarstedt S-Monovette® (15 x 75 mm), (4; 4,3; 5,5 ml)
[28]	*	Sarstedt S-Monovette® (15 x 92 mm), (7,5; 8,2; 8,5 ml)
[28]	*	10 ml Thermo Nalgene® Oak Ridge (16 x 81,5 mm)
[28]	15046	14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt® 14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®
[28]	15118	10 ml probówka szklana (16 x 100 mm) 10 ml glass tube (16 x 100 mm)

14020

[20]	15053	10 ml probówka z pokrywką (16 x 106 mm) 10 ml tube with cap (16 x 106 mm)
[20]	*	Sarstedt S-Monovette® (15 x 75 mm), (4; 4,3; 5,5 ml)
[20]	*	10 ml Thermo Nalgene® Oak Ridge (16 x 81,5 mm)
[20]	15046	14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt® 14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®
[20]	15121	10 ml probówka z dnem okrągłym i pokrywką (17 x 70 mm) 10 ml tube, round bottom, with cap (17 x 70 mm)
[20]	*	BD Vacutainer® (16 x 100 mm), (2,5-11 ml)
[20]	*	Greiner Vacuette® (16 x 100 mm), (7-9 ml)
[20]	*	Sarstedt S-Monovette® (15 x 92 mm), (7,5; 8,2; 8,5 ml)
[20]	*	Sarstedt S-Monovette® (16 x 92 mm), (9; 10 ml)
[20]	15118	10 ml probówka szklana (16 x 100 mm) 10 ml glass tube (16 x 100 mm)
[20]	*	13 ml probówka (16x100mm), Sarstedt® nr 62.515.006 13 ml tube (16 x 100 mm), Sarstedt® no. 62.515.006

14021

[40]	*	2-1,5 ml probówka (10,8x41,8 mm), Eppendorf®; [15011], 2 ml (10,8x41,8 mm); [15128], 1,5ml (10,8x40,5 mm)
[40]	*	2 ml probówki z filtrem - spin columns (10,8 x 46 mm) 2 ml spin columns (with filter) (10,8 x 46 mm); [15011], 2 ml (10,8x41,8 mm); [15128], 1,5ml

14023

[4]	15055	30 ml probówka z pokrywką (25,4 x 103,2 mm) 30 ml tube with cap (25,4 x 103,2 mm)
[4]	15222	30 ml probówka z pokrywką (25 x 94mm), Sterilin® 30 ml tube with cap (25 x 94 mm), Sterilin®
[4]	15223	30 ml probówka z pokrywką (25 x 94 mm), Sterilin® 30 ml tube with cap (25 x 94 mm), Sterilin®
[4]	*	28 ml Thermo Nalgene® Oak Ridge (25,4 x 101,8 mm)

* probówka niedostępna w ofercie MPW lub dostępny odpowiednik (np.[15050]), patrz kolumna z prawej
tube is not offered by MPW or equivalent is available (e.g. [15050]), see column on the right

A. Wyposażenie dodatkowe/Optional accessories
MPW-260/R/RH

[4]	*	30 ml probówka z pokrywką (25,5 x 94 mm), Nalgene® 30 ml tube with cap (25,5 x 94 mm), Nalgene®
[4]	*	30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm) 30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm)
[4]	*	25 ml probówka szklana (25 x 100 mm) 25 ml glass tube (25 x 100 mm)
14026		
[4]	*	50 ml probówka z dnem stożkowym z rantem (30 x 115 mm), Greiner® 50 ml tube, conical bottom, skirted (30 x 115 mm), Greiner®
14026+14188		
[4]	15051	50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm) 50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm)
[4]	*	50 ml probówka z dnem stożkowym z zakrętką (30 x 117 mm), Falcon®; [15052] 50ml (30 x 117mm) 50 ml tube, conical bottom, with cap (30 x 117 mm), Falcon®; [15052] 50ml Sarstedt® (30 x 117 mm)
[4]	*	50 ml probówka z dnem stożkowym bez rantu (30 x 115 mm), Greiner® 50 ml tube, conical bottom, without skirt (30 x 115 mm), Greiner®
[4]	*	50 ml probówka Advanced Oak Ridge (29x102 mm), Herolab® nr 25 32 11 50 ml tube, Advanced Oak Ridge (29 x 102 mm), Herolab® no. 25 32 11
14028		
[4]	*	50 ml probówka szklana (35 x 100 mm) 50 ml glass tube (35 x 100 mm)
14029		
[48]	*	Sarstedt S-Monovette® (11 x 92 mm), (4,5; 5 ml)
[48]	*	Sarstedt S-Monovette® (11 x 66 mm), (1,6; 2; 2,7; 3; 3,1 ml)
[48]	*	7 ml probówka szklana (12 x 100 mm) 7 ml glass tube (12 x 100 mm)
[48]	*	5 ml probówka szklana (12 x 75 mm) 5 ml glass tube (12 x 75 mm)
[48]	*	5 ml probówka z korkiem (12 x 85 mm), Sarstedt® 5 ml tube with cap (12 x 85 mm), Sarstedt®
[48]	*	6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt® 6 ml tube with cap (11,5 x 92 mm), Sarstedt®
14100+14196		
[4]	15040	100 ml probówka z pokrywką (45,2 x 103,7 mm) 100 ml tube with cap (45,2 x 103,7 mm)
14027		
[4]	*	15 ml probówka z dnem stożkowym z zakrętką (17 x 120 mm), Falcon®; [15050], 15ml (17 x 120 mm) 15 ml tube, conical bottom, with cap (17 x 120 mm), Falcon®; [15050] 15ml Sarstedt®(17 x 120 mm)
14100+14188		
[4]	*	100 ml probówka szklana (44 x 100 mm) 100 ml glass tube (44 x 100 mm)
13201+17202		
14013		
[32]	*	BD Vacutainer® (13 x 75 mm), (1,6-5,3 ml)
[32]	*	BD Vacutainer® (13 x 100 mm), (4-7 ml)
[32]	*	Greiner Vacuette® (13 x 75 mm), (1-4,5 ml)
[32]	*	Greiner Vacuette® (13 x 100 mm), (3,5-6 ml)
[32]	*	Sarstedt S-Monovette® (13 x 65 mm), (2,6; 2,9; 3,4; 3,8 ml)
[32]	*	Sarstedt S-Monovette® (13 x 75 mm), (2,7; 3; 4,3 ml)
[32]	*	Sarstedt S-Monovette® (11 x 92 mm), (4,5; 5 ml)
[32]	*	Sarstedt S-Monovette® (13 x 90 mm), (4,9; 5,6 ml)
[32]	*	7 ml probówka szklana (12 x 100 mm) 7 ml glass tube (12 x 100 mm)
[32]	*	5 ml probówka szklana (12 x 75 mm) 5 ml glass tube (12 x 75 mm)
[32]	*	5 ml probówka z korkiem (12 x 85 mm), Sarstedt® 5 ml tube with cap (12 x 85 mm), Sarstedt®
[32]	*	6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt® 6 ml tube with cap (11,5 x 92 mm), Sarstedt®
14016		
[28]	15053	10 ml probówka z pokrywką (16 x 106 mm) 10 ml tube with cap (16 x 106 mm)
[28]	*	BD Vacutainer® (16 x 100 mm), (2,5-11 ml)
[28]	*	Greiner Vacuette® (16 x 100 mm), (7-9 ml)
[28]	*	Sarstedt S-Monovette® (15 x 75 mm), (4; 4,3; 5,5 ml)
[28]	*	10 ml Thermo Nalgene® Oak Ridge (16 x 81,5 mm)
[28]	15046	14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt® 14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®
[28]	15118	10 ml probówka szklana (16 x 100 mm) 10 ml glass tube (16 x 100 mm)
14020		
[20]	15053	10 ml probówka z pokrywką (16 x 106 mm) 10 ml tube with cap (16 x 106 mm)
[20]	*	Sarstedt S-Monovette® (15 x 75 mm), (4; 4,3; 5,5 ml)
[20]	*	10 ml Thermo Nalgene® Oak Ridge (16 x 81,5 mm)

* probówka niedostępna w ofercie MPW lub dostępny odpowiednik (np:[15050]), patrz kolumna z prawej
tube is not offered by MPW or equivalent is available (e.g. [15050]), see column on the right

A. Wyposażenie dodatkowe/Optional accessories
MPW-260/R/RH

[20]	15046	14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt® 14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®
[20]	15121	10 ml probówka z dnem okrągłym i pokrywką (17 x 70 mm) 10 ml tube, round bottom, with cap (17 x 70 mm)
[20]	*	BD Vacutainer® (16 x 100 mm), (2,5-11 ml)
[20]	*	Greiner Vacuette® (16 x 100 mm), (7-9 ml)
[20]	15118	10 ml probówka szklana (16 x 100 mm) 10 ml glass tube (16 x 100 mm)
[20]	*	13 ml probówka (16x100mm), Sarstedt® nr 62.515.006 13 ml tube (16 x 100 mm), Sarstedt® no. 62.515.006
14021		
[40]	*	2-1,5 ml probówka (10,8x41,8 mm), Eppendorf®; [15011], 2 ml (10,8x41,8 mm); [15128], 1,5ml (10,8x40,5 mm)
[40]	*	2 ml probówki z filtrem - spin columns (10,8 x 46 mm) 2 ml spin columns (with filter) (10,8 x 46 mm); [15011], 2 ml (10,8x41,8 mm); [15128], 1,5ml
14023		
[4]	15055	30 ml probówka z pokrywką (25,4 x 103,2 mm) 30 ml tube with cap (25,4 x 103,2 mm)
[4]	15222	30 ml probówka z pokrywką (25 x 94mm), Sterilin® 30 ml tube with cap (25 x 94 mm), Sterilin®
[4]	15223	30 ml probówka z pokrywką (25 x 94 mm), Sterilin® 30 ml tube with cap (25 x 94 mm), Sterilin®
[4]	*	28 ml Thermo Nalgene® Oak Ridge (25,4 x 101,8 mm)
[4]	*	30 ml probówka z pokrywką (25,5 x 94 mm), Nalgene® 30 ml tube with cap (25,5 x 94 mm), Nalgene®
[4]	*	30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm) 30 ml Thermo Nalgene® Oak Ridge (25,5 x 94,3 mm)
[4]	*	25 ml probówka szklana (25 x 100 mm) 25 ml glass tube (25 x 100 mm)
14026+14188		
[4]	15051	50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm) 50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm)
14028		
[4]	*	50 ml probówka szklana (35 x 100 mm) 50 ml glass tube (35 x 100 mm)
14029		
[48]	*	Sarstedt S-Monovette® (11 x 66 mm), (1,6; 2; 2,7; 3; 3,1 ml)
[48]	*	7 ml probówka szklana (12 x 100 mm) 7 ml glass tube (12 x 100 mm)
[48]	*	5 ml probówka szklana (12 x 75 mm) 5 ml glass tube (12 x 75 mm)
[48]	*	5 ml probówka z korkiem (12 x 85 mm), Sarstedt® 5 ml tube with cap (12 x 85 mm), Sarstedt®
[48]	*	6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt® 6 ml tube with cap (11,5 x 92 mm), Sarstedt®
14100+14196		
[4]	15040	100 ml probówka z pokrywką (45,2 x 103,7 mm) 100 ml tube with cap (45,2 x 103,7 mm)
14100+14188		
[4]	*	100 ml probówka szklana (44 x 100 mm) 100 ml glass tube (44 x 100 mm)
13201+17203		
14021		
[40]	*	2-1,5 ml probówka (10,8x41,8 mm), Eppendorf®; [15011], 2 ml (10,8x41,8 mm); [15128], 1,5ml (10,8x40,5 mm)
[40]	*	2 ml probówki z filtrem - spin columns (10,8 x 46 mm) 2 ml spin columns (with filter) (10,8 x 46 mm); [15011], 2 ml (10,8x41,8 mm); [15128], 1,5ml
14026		
[4]	*	50 ml probówka z dnem stożkowym z rantem (30 x 115 mm), Greiner® 50 ml tube, conical bottom, skirted (30 x 115 mm), Greiner®
14026+14188		
[4]	15051	50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm) 50 ml Thermo Nalgene® Oak Ridge (28,8 x 106,7 mm)
[4]	*	50 ml probówka z dnem stożkowym z zakrętką (30 x 117 mm), Falcon®; [15052] 50ml (30 x 117mm) 50 ml tube, conical bottom, with cap (30 x 117 mm), Falcon®; [15052] 50ml Sarstedt® (30 x 117 mm)
[4]	*	50 ml probówka z dnem stożkowym bez rantu (30 x 115 mm), Greiner® 50 ml tube, conical bottom, without skirt (30 x 115 mm), Greiner®
[4]	*	50 ml probówka Advanced Oak Ridge (29x102 mm), Herolab® nr 25 32 11 50 ml tube, Advanced Oak Ridge (29 x 102 mm), Herolab® no. 25 32 11
14028		
[4]	*	50 ml probówka szklana (35 x 100 mm) 50 ml glass tube (35 x 100 mm)
14100+14196		
[4]	15040	100 ml probówka z pokrywką (45,2 x 103,7 mm) 100 ml tube with cap (45,2 x 103,7 mm)
14027		

* probówka niedostępna w ofercie MPW lub dostępny odpowiednik (np.[15050]), patrz kolumna z prawej
tube is not offered by MPW or equivalent is available (e.g. [15050]), see column on the right

A. Wyposażenie dodatkowe/Optional accessories
MPW-260/R/RH

[4] * 15 ml probówka z dnem stożkowym z zakrętką (17 x 120 mm), Falcon®; [15050], 15ml (17 x 120 mm)
15 ml tube, conical bottom, with cap (17 x 120 mm), Falcon®; [15050] 15ml Sarstedt®(17 x 120 mm)
14100+14188

[4] * 100 ml probówka szklana (44 x 100 mm)
100 ml glass tube (44 x 100 mm)

13215
14815 Rmax 138 RCF 2469

[8] 15121 10 ml probówka z dnem okrągłym i pokrywką (17 x 70 mm)
10 ml tube, round bottom, with cap (17 x 70 mm)

[8] * Sarstedt S-Monovette® (15 x 75 mm), (4; 4,3; 5,5 ml)

[8] * 10 ml Thermo Nalgene® Oak Ridge (16 x 81,5 mm)

14082+14815 Rmax 138 RCF 2469

[8] * BD Vacutainer® (13 x 75 mm), (1,6-5,3 ml)

[8] * Greiner Vacuette® (13 x 75 mm), (1-4,5 ml)

[8] * Sarstedt S-Monovette® (11 x 66 mm), (1,6; 2; 2,7; 3; 3,1 ml)

[8] * Sarstedt S-Monovette® (13 x 65 mm), (2,6; 2,9; 3,4; 3,8 ml)

[8] * Sarstedt S-Monovette® (13 x 75 mm), (2,7; 3; 4,3 ml)

[8] * Sarstedt V-Monovette urine tube (13 x 75 mm)

[8] * BD urine tube (13 x 75 mm)

[8] * 5 ml probówka szklana (12 x 75 mm)

5 ml glass tube (12 x 75 mm)

[8] * 5 ml probówka z korkiem (12 x 85 mm), Sarstedt®

5 ml tube with cap (12 x 85 mm), Sarstedt®

[8] * 6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt®

6 ml tube with cap (11,5 x 92 mm), Sarstedt®

13113 R max 121 RCF 2164
bez wkładki/without adapter Rmax 121 RCF 2164

[48] * BD Vacutainer® (13 x 75 mm), (1,6-5,3 ml)

[48] * Greiner Vacuette® (13 x 75 mm), (1-4,5 ml)

[48] * Sarstedt S-Monovette® (11 x 66 mm), (1,6; 2; 2,7; 3; 3,1 ml)

[48] * Sarstedt S-Monovette® (13 x 65 mm), (2,6; 2,9; 3,4; 3,8 ml)

13215 R max 138 RCF 2469
bez wkładki/without adapter Rmax 138 RCF 2469

[8] 15046 14 ml probówka z pokrywką (16,8 x 113,7 mm), Sarstedt®

14 ml tube with cap (16,8 x 113,7 mm), Sarstedt®

[8] 15053 10 ml probówka z pokrywką (16 x 106 mm)

10 ml tube with cap (16 x 106 mm)

[8] * BD Vacutainer® (16 x 100 mm), (2,5-11 ml)

[8] * Greiner Vacuette® (16 x 100 mm), (7-9 ml)

[8] * Sarstedt S-Monovette® (15 x 92 mm), (7,5; 8,2; 8,5 ml)

[8] * Sarstedt S-Monovette® (16 x 92 mm), (9; 10 ml)

[8] * 15 ml probówka z dnem stożkowym z zakrętką (17 x 120 mm), Falcon®; [15050], 15ml (17 x 120 mm)

15 ml tube, conical bottom, with cap (17 x 120 mm), Falcon®; [15050] 15ml Sarstedt®(17 x 120 mm)

[8] 15118 10 ml probówka szklana (16 x 100 mm)

10 ml glass tube (16 x 100 mm)

[8] * 15 ml Thermo Nalgene® (16 x 113 mm)

15 ml Thermo Nalgene® (16 x 113 mm)

14082 Rmax 138 RCF 2469

[8] * BD Vacutainer® (13 x 100 mm), (4-7 ml)

[8] * Greiner Vacuette® (13 x 100 mm), (3,5-6 ml)

[8] * Sarstedt S-Monovette® (11 x 92 mm), (4,5; 5 ml)

[8] * 7 ml probówka szklana (12 x 100 mm)

7 ml glass tube (12 x 100 mm)

[8] * 6 ml probówka z pokrywką (11,5 x 92 mm), Sarstedt®

6 ml tube with cap (11,5 x 92 mm), Sarstedt®

12218

RPM 3000 RCF 916 Rmax 91 \neq 90

13219
bez wkładki/without adapter

[2] * płytki titracyjna MTP 28,8ml (86x128x15/17,5 mm)

microtiter plate MTP 28,8 ml (86 x 128 x 15/17,5 mm)

12300

RPM 13000 RCF 16816 Rmax 89 \neq 90

bez pojemnika/without bucket
bez wkładki/without adapter

[24] * 37 μ l kapilara hematokrytowa (1,4 x 75 mm)

37 μ l micro-hematocrit capillary tube (1,4 x 75 mm)

12300C

RPM 13000 RCF 16816 Rmax 89 \neq 90

bez pojemnika/without bucket
bez wkładki/without adapter

* probówka niedostępna w ofercie MPW lub dostępny odpowiednik (np:[15050]), patrz kolumna z prawej
tube is not offered by MPW or equivalent is available (e.g. [15050]), see column on the right

A. Wyposażenie dodatkowe/Optional accessories**MPW-260/R/RH**

[24] * 37 µl kapilara hematokrytowa (1,4 x 75 mm)
37 µl micro-hematocrit capillary tube (1,4 x 75 mm)

Suma końcowa

EU DECLARATION OF CONFORMITY

This EU declaration of conformity is issued under the sole responsibility of the manufacturer.

Manufacturer: "MPW MED. INSTRUMENTS" SPÓŁDZIELNIA PRACY
46 Boremlowska Street, 04-347 Warsaw, Poland

The Quality Management System complies with the standards: PN-EN ISO 9001:2015, PN-EN ISO 13485:2016

SRN: PL-MF-000032831

Device name: Laboratory centrifuge MPW-260
(with the accessory indicated in the operating instructions provided with the centrifuge)

BASIC UDI-DI: 590538636-IVD-CEN-005-69

Catalogue numbers: 10260/2-56 10260/1-56 10260/1-56/100
10260/1-56/110 10260/1-56/127

The aforementioned device is in conformity with the following EU regulations and directives:

2017/746 (IVDR) REGULATION (EU) 2017/746 OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL of 5 April 2017 on in vitro diagnostic medical devices and repealing Directive 98/79/EC and Commission Decision 2010/227/EU, including the changes published prior to the date of this declaration.

2011/65/EU (RoHS 2) DIRECTIVE 2011/65/EU OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, including the changes published prior to the date of this declaration.

Intended purpose: The device is intended for the separation of the mixtures of the liquid substances derived from the human body, including blood, urine, and other body fluids, and for the preparation of the samples intended for further in vitro diagnostics procedures.

Risk class: Class A
(in accordance with the rule 5 of Annex VIII of Regulation (EU) 2017/746).

The conformity assessment of the device and accessory has been carried out in accordance with Article 48(10) of Regulation (EU) 2017/746.

Wojciech Anisiewicz
Vice-President of the Management Board

Łukasz Szański
President of the Management Board

EU DECLARATION OF CONFORMITY

This EU declaration of conformity is issued under the sole responsibility of the manufacturer.

Manufacturer: "MPW MED. INSTRUMENTS" SPÓŁDZIELNIA PRACY
46 Boremlowska Street, 04-347 Warsaw, Poland

The Quality Management System complies with the standards: PN-EN ISO 9001:2015, PN-EN ISO 13485:2016

SRN: PL-MF-000032831

Device name: Refrigerated laboratory centrifuge MPW-260R
(with the accessory indicated in the operating instructions provided with the centrifuge)

BASIC UDI-DI: 590538636-IVD-CEN-006-6C

Catalogue numbers: 10260R/2-5 10260R/2-6 10260R/1-6
10260R/1-6/100 10260R/1-6/110 10260R/1-6/127

The aforementioned device is in conformity with the following EU regulations and directives:

2017/746 (IVDR) REGULATION (EU) 2017/746 OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL of 5 April 2017 on in vitro diagnostic medical devices and repealing Directive 98/79/EC and Commission Decision 2010/227/EU, including the changes published prior to the date of this declaration.

2011/65/EU (RoHS 2) DIRECTIVE 2011/65/EU OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, including the changes published prior to the date of this declaration.

Intended purpose: The device is intended for the separation of the mixtures of the liquid substances derived from the human body, including blood, urine, and other body fluids, and for the preparation of the samples intended for further in vitro diagnostics procedures.

Risk class: Class A
(in accordance with the rule 5 of Annex VIII of Regulation (EU) 2017/746).

The conformity assessment of the device and accessory has been carried out in accordance with Article 48(10) of Regulation (EU) 2017/746.

Wojciech Anisiewicz
Vice-President of the Management Board


Lukasz Salański
President of the Management Board

DECLARATION OF DECONTAMINATION

(repair)

In order to protect our employees please fill out the declaration of decontamination completely before sending centrifuge to the manufacturer (repair).

If it is impossible to completely and effectively decontaminate the device, it should be treated in accordance with the regulations for medical waste.

1. Device:

– type:

– serial No.:

2. Description of decontamination

(see user manual)

.....

.....

.....

.....

3. Decontamination carried out by:

name:

4. Date and signature:

.....

DECLARATION OF DECONTAMINATION

(return)

In order to protect our employees please fill out the declaration of decontamination completely before sending centrifuge to the manufacturer (return).

If it is impossible to completely and effectively decontaminate the device, it should be treated in accordance with the regulations for medical waste.

1. Device:

– type:

– serial No.:

2. Description of decontamination

(see user manual)

.....

.....

.....

.....

3. Decontamination carried out by:

name:

4. Date and signature:

.....

NOMOGRAM

Centrifuging radius [cm]

R.C.F. (x "g")
multiple of
gravitational
acceleration

[r.p.m.]

